





# TR 9650

## Mobile Bath Lift System

Hygiene Equipment  
Designed for People

### Quick guide for the Chair Lift

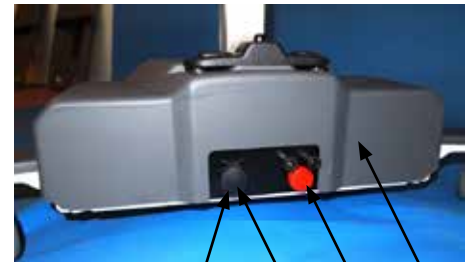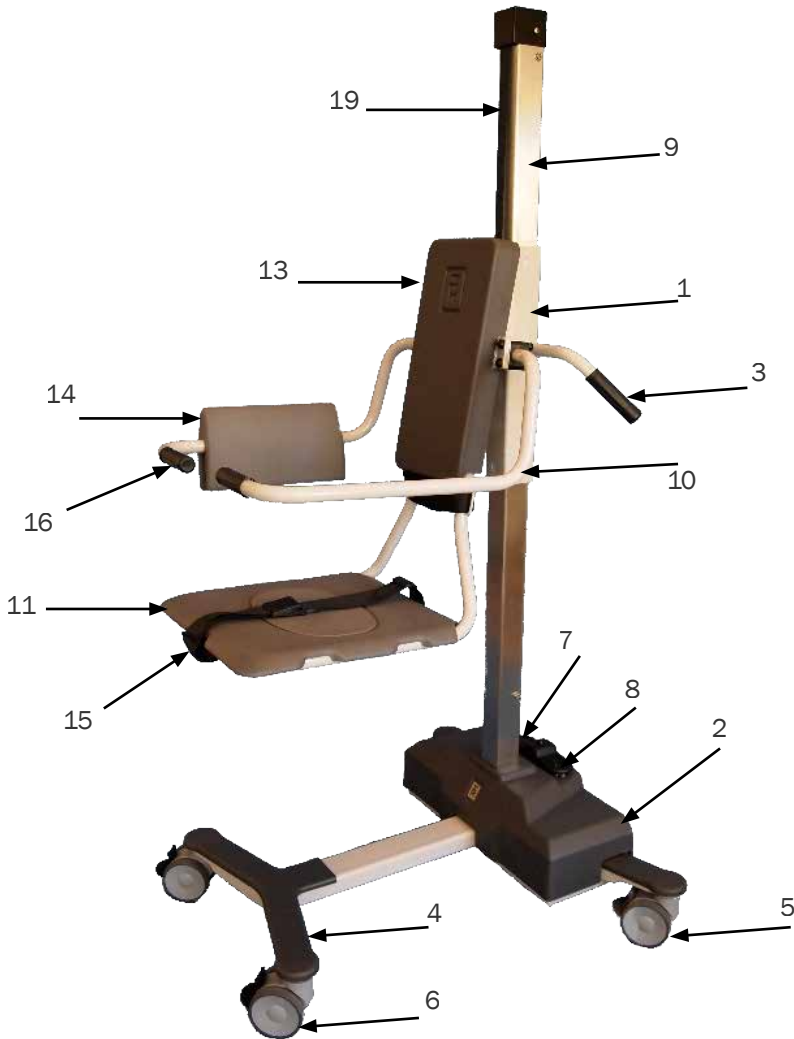

Open lid



20  
21 22

**Working Load Limit: 150 kg (330 lbs)**

Position	Description		
		10	Rotating arm/back supports
1	Support column	11	Hygiene seat and slot
2	Chassis Base	13	"TR" Protection pad
3	Transport handle	14	Backrest Pad
4	Chassis leg	15	Safety belts
5	Positive brake castors, 100mm Twin	16	Grip
6	Straight steering castors, 100mm Twin	19	Drive belt
7	Foot control for raising	20	Low battery audible tone
8	Foot control for lowering	21	Battery charger input
9	Stainless Steel column	22	Emergency Stop



**TR EQUIPMENT**



# TR 9650

## Mobile Bath Lift System

Hygiene Equipment  
Designed for People



### Quick guide for the Chairlift

- Only qualified staff, who understand all operating and safety procedures may operate the equipment.
- Safety belts must be securely fastened around patient at all times.
- Never leave the patient unattended.
- Castor brakes must be securely locked on the Lifter before transferring a patient on to the chair.
- The patient must be positioned comfortably on the chair, with folded rotating arm.
- Make sure that all parts that come in contact with the patient are intact and without sharp edges.
- Ensure that the Lift can be moved, raised or lowered safely to avoid any pinch injuries to patient.
- Operate the Lift carefully when transporting patients through narrow corridors or over irregular uneven surfaces.
- The TR battery charger is the only one approved for charging.
- Do not use battery charger in wet or damp areas - potential shock hazard.
- It is recommended that the batteries are charged daily or whenever Lifter is not in use. Batteries can not be overcharged.
- Check that castor brakes and straight steering are free from damage.
- Be sure that all details are in place and free from damage and all screws are in place.
- Pads must be kept well maintained in place.
- Lift and chair must be inspected visually for sign of wear. Lift should not be used if defects are noted. Prompt service is required.

**Warning:** Never use the Lifter in any way that varies from the safety precautions and operating procedures detailed in the operating manual. The equipment must not be operated by anyone until the user has fully read and understands the instructions in the operating manual.

#### **Safety autostop mechanism:**

If the chair hits an obstacle when lowering, the Lift will safely stop automatically. To release it, just press the "up" right arrow foot control. Same procedure if the lifting belt is slacking. If lifting belt fails, the chair will stop/lock automatically on the column. Call for technical assistance.

**For further information see Operating manual**



TR EQUIPMENT