




INVACARE®
RELIANT™ 350
 STAND-UP LIFT



Yes, you can.®

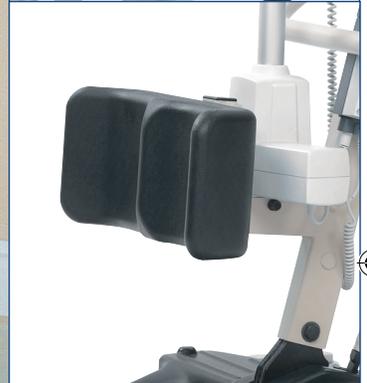

The **Invacare® Reliant™ Series Stand-Up Lifts** are an ideal product for everyday patient handling, ensuring caregiver security and patient peace of mind during the transfer. Intended for use with weight-bearing patients as well as those needing rehabilitation support, the Invacare Reliant Stand-Up lift will be a valuable addition to your safety management system.



Features

Innovative adjustment features adapt to a wide variety of body sizes and shapes

Low friction casters significantly improve rollability and contribute to caregiver safety

Multi-function slings allow for quick toileting, stand assist or full-support seated transfer

Adjustable leg support features a soft pad with lateral contour that secures legs in three positions for comfort at any height

Wide, non-slip footplate provides stable base for transfer

Under-bed clearance of 4.5" allows accessibility to most nursing home, hospital and homecare bed styles

Available with power opening base (model # RPS350-2)

Removable battery eliminates downtime


Ergonomic features make it ideal for toilet transfers, as well as for dependent and rehabilitation purposes

Features include fast, simple sling attachment that protects against accidental disengagement; wide push handle; good access for clothing removal; and comfortable slings with extended back support that prevent pressure in the under-arm area

Adjustable leg strap option secures leg against knee pad for optimum positioning on footplate

Manufactured in an ISO13485 facility

HCPCS Code: E0635

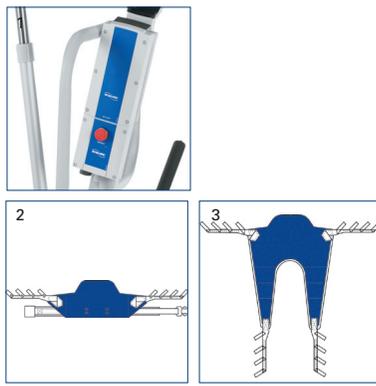

1. Electronics
2. Standing sling* Models # R130 and R131
3. Transfer sling* Models # R134 and R136

*Complete instructions for lift usage and sling procedures are found in the Invacare operating manual. Only Invacare® slings and accessories should be used on Invacare® lifts.



Invacare Corporation
www.invacare.com

USA
One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada
570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

Model #	Height At Sling Hook Up	Base Dimensions	Casters	Weight	Weight Capacity	Low Battery Alarm	Emergency Controls
RPS350-1	39.6"-63.7" H	4.5"H Clearance 36.8"W Open 25.7"W Closed 35.5"L	3" Front 5" Rear	108 lb.	350 lb.	Audio	Stop and Down
RPS350-2	39.6"-63.7" H	4.5"H Clearance 36.8"W Open 25.7"W Closed 35.5"L	3" Front 5" Rear	108 lb.	350 lb.	Audio	Stop and Down

	Battery	Charge Time	Charger Input	Charger Output	Limited Warranty
RPS350-1	24V DC	6 hr. Max	100 - 240V AC	29.5V DC	3 Years Lift 1 Year Electronics
RPS350-2	24V DC	6 hr. Max	100 - 240V AC	29.5V DC	3 Years Lift 1 Year Electronics

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.

