



CS5 bed featured Belmont Bed Ends

CARROLL CS SERIES™
CS5™ BED



Yes, you can.®

It's time to change your definition of value. The **Carroll CS Series™ CS5™ Bed** was developed for facilities requiring mobility at any height at a great value, and who appreciate the trademark Carroll CS Series Bed design elegance including the slat deck and tool-less assembly.

Invacare® beds, mattresses and rails, excluding therapeutic support surfaces and bariatric beds, when used as a system and sold new and unused after August 1, 2007, meet the recommended FDA Guidelines for bed rail entrapment. Visit www.invacare.com/bedsafety for additional safety and entrapment information.

Features

Mobility at any Height – The CS5 bed is mobile in any position providing convenience for the caregiver.

Pure Vertical Rise – Raises straight up and down throughout the travel range, with zero horizontal movement of the casters minimizing potential damage to the wall or flooring.

Exclusive Slat Deck Design – Provides a comfortable sleep surface, allows for improved infection control and helps reduce wear and tear on the mattress.


Ergonomic Hand Pendant – Water-resistant¹ and rugged hand pendant designed to fit comfortably in the resident's hand.

Dual Pendant Connections – Pendant connects to either side of the bed, providing customized convenience for the resident. Pendant cord is attached to frame, and will not interfere with lifts or tables.

Tool-less Assembly – The CS5 bed features tool-less assembly, minimizing set-up time.

Customizable Options –

- One-touch floor lock
- 39" or 42" Wide Deck
- 4" Bed Extender (84" Length)
- Attendant Control Panel
- Steering Lock
- Hi-Lo Lockout
- On board bathing
- 4 rail styles

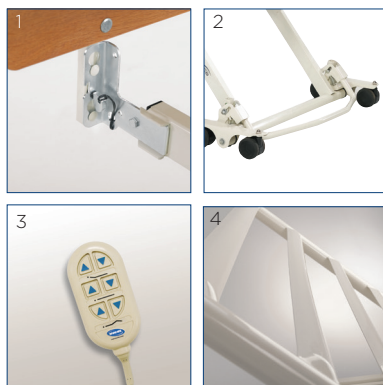

1. Tool-less Assembly 3. Ergonomic Hand Pendant
2. Mobility at any Height 4. Slat Deck Design



Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

©2020

Invacare Corporation.
All rights reserved.
Trademarks are identified by the symbols TM and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.
Form No. 20-031
200203

Model #	Height Range	Sleep Surface	Weight Capacity*	Regulation Compliance	Limited Warranty
IHCS5	9" - 30" H	76" - 80" L 36", 39", 42" W	500 lb.	†Classified to: UL 60601-1 IEC 60601-2-38 CAN/CSA C22.2 No. 601.1	10 Years Frame 5 Years Mechanical 3 Years Electronics

NOTE: Whenever bed rails are used as part of any hospital or homecare bed system, the potential for bed rail entrapment exists. To reduce this risk, a patient must always be assessed regarding their risk of entrapment before any product is used. Routine patient monitoring is also required once a product is in use. To learn more about the risk of bed rail entrapment please visit the US FDA's website at <http://www.fda.gov> and search for "Bed Rail Entrapment."

* Maximum safe working load.

1. Pendant is water-resistant to IPX6 per IEC 60529.

†Testing conducted to select portions of UL/IEC Standards.

