



Molift

Hoists, slings & transfer aids
Catalogue

 **etac**®

A photograph of a caregiver, a young woman with brown hair in a ponytail wearing a blue shirt, assisting an elderly woman. The elderly woman has short brown hair, wears glasses, an orange sweater, and blue jeans. She is holding a grey walker. They are both smiling and looking at each other in a bright, indoor setting with a window in the background.

Etac - Creating Possibilities

Etac is a world-leading developer and provider of ergonomic assistive devices and patient handling equipment. Through our specialized product brands, we offer state-of-the-art products for a wide range of daily needs and care settings, for all stages of life. Whether the solution is for a child or an elderly person, our aim is to promote the individual's abilities and improve the caregiver's working conditions – through sustainable quality, functionality and design.

molift®

Molift by Etac - offers ergonomic hoists, slings and transfer assistive devices ensuring safe and comfortable transfers in various care settings with both the caregiver and patient in mind.



Sit-to-stand solutions

Molift Raiser Pro.....	7
Etac Turner Pro.....	8
Molift Quick Raiser 205.....	9
Molift RgoSling Active	10



Mobile Hoists

Molift Smart 150.....	12
Molift Partner 255	13
Molift Mover	14
Sling bars for Molift mobile hoists	16
Molift Service Tool.....	17



Ceiling Hoists

Molift Nomad.....	19
Molift Air	20
Molift Rail System	22
Sling bars for Molift ceiling hoists.....	24
Molift Quattro	25
Molift Duo	26

Stretchers



Molift RgoSling Fabric Stretcher	28
Molift Stretcher	29

Molift RgoSling 4-point



Molift RgoSling HighBack	32
Molift RgoSling MediumBack.....	32
Molift RgoSling Toilet	33
Molift RgoSling Amputee	33
Molift RgoSling Shadow.....	34
Molift RgoSling Comfort HighBack.....	34
Molift RgoSling Ambulating Vest.....	35
Molift RgoSling Fabric Stretcher	35
Molift RgoSling MediumBack Plus.....	37
Molift RgoSling Repositioning sheet.....	37

Molift EvoSling 2-point



Molift EvoSling HighBack.....	40
Molift EvoSling MediumBack	40
Molift EvoSling LowBack.....	41
Molift EvoSling Hygiene	41
Molift EvoSling Ampu MediumBack.....	42
Molift EvoSling Comfort MediumBack.....	43
Molift EvoSling FlexiStrap	43

Molift UnoSling

Molift UnoSling HighBack.....	45
Facts.....	46
Molift Slings accessories.....	47

Molift products are designed for life


Molift products are designed to give users a better quality of life while providing caregivers with the best support and optimum work environment. Our objective is to make solutions that improve the quality of life each day for users and give caregivers a safer working environment, enabling them to focus more of their time and energy on the users.

Best in class - Lasting durability

Molift products give disabled people the opportunity to make the most of life - go to school, work and travel – a better life through safe, efficient and comfortable transfer solutions. Molift transfer solutions are designed to be the best in class when it comes to product quality, reliability and durability.

The Molift story

Molift became a part of Etac 2010. The brand was founded more than 30 years ago in Norway by Hans Kasper Andresen. A Norwegian inventor, engineer and designer who, through his own experience (he was hospitalised due to a heart attack), saw the need for lifting and moving users in a secure way.

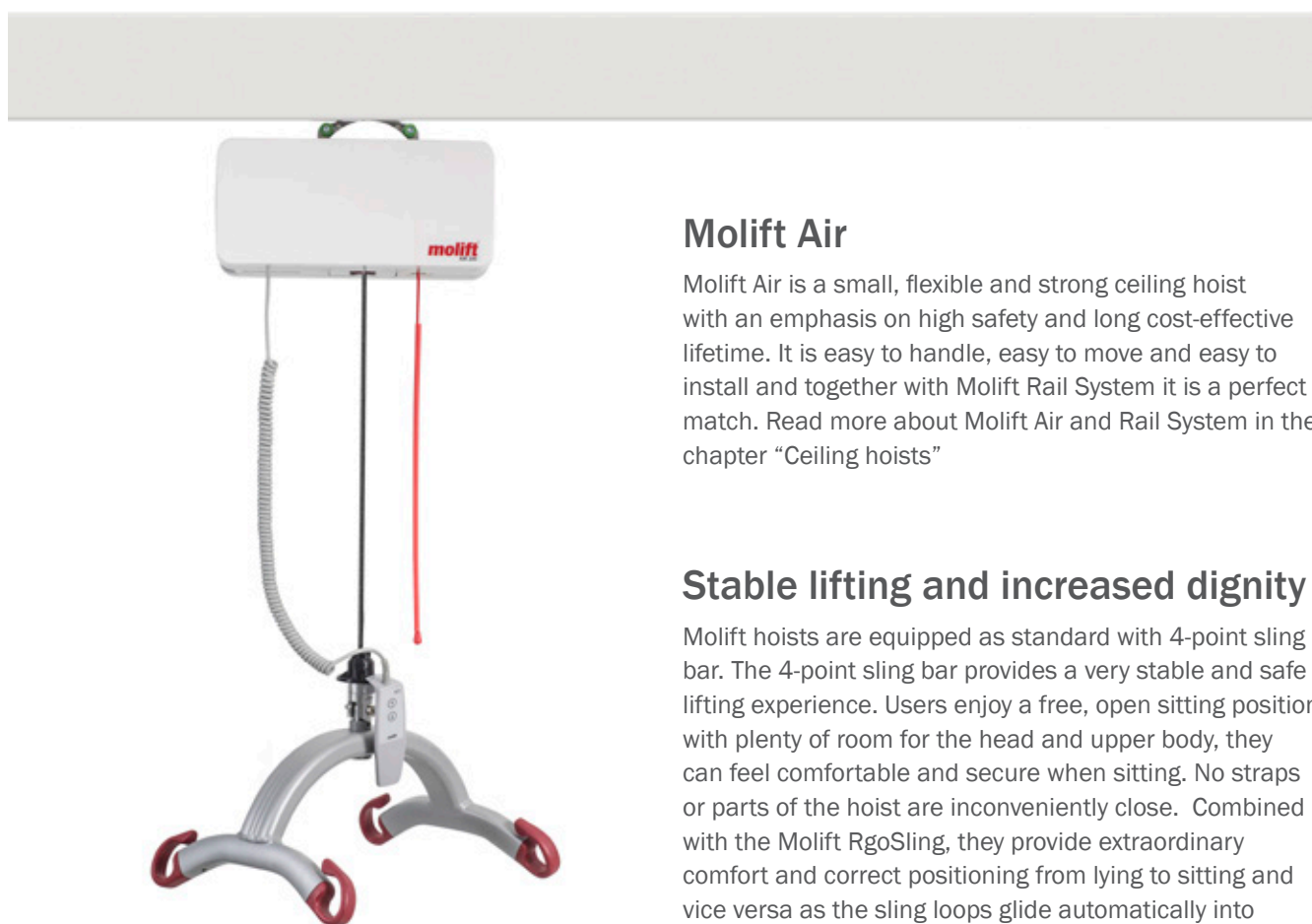

Molift Air

Molift Air is a small, flexible and strong ceiling hoist with an emphasis on high safety and long cost-effective lifetime. It is easy to handle, easy to move and easy to install and together with Molift Rail System it is a perfect match. Read more about Molift Air and Rail System in the chapter “Ceiling hoists”

Stable lifting and increased dignity

Molift hoists are equipped as standard with 4-point sling bar. The 4-point sling bar provides a very stable and safe lifting experience. Users enjoy a free, open sitting position with plenty of room for the head and upper body, they can feel comfortable and secure when sitting. No straps or parts of the hoist are inconveniently close. Combined with the Molift RgoSling, they provide extraordinary comfort and correct positioning from lying to sitting and vice versa as the sling loops glide automatically into position during hoisting. Most of the Molift hoists are also equally available with 2-point sling bars.



Secure, efficient and outstanding

With the increasing need for efficiency in the healthcare sector and a growing focus on users' and caregivers' needs and safety, top class transfer products are essential.

Cost effective solutions

Studies show that manual lifting and other user-handling tasks are high-risk activities for both carers and users. The aim for Molift products is to give the best and most secure user care, avoid heavy lifting situations which can put caregivers' backs at risk, while at the same time as providing the healthcare sector with the most cost effective solutions. Research and studies over for the past 20 years have shown that the costs of proper equipment are far lower than those associated with carers' work related injuries.

Our standards

Molift is certified according to the ISO standards for our trade. This affects how we drive and safeguard the quality of our products and our processes.

Systems for quality assurance

We are certified according to ISO 9001, which means that all our work processes, from product development and production to the products and their delivery meets the most stringent quality requirements.

Traceability

Molift meets all the requirements of ISO 13485 with regard to traceability of our products. All our hoists and slings have bar codes and unique serial numbers to facilitate traceability for our distributors on their market.

Environmental requirements

We work continuously to minimise the impact of our operation on the environment, and we are certified according to ISO 14001. This standard contains legal requirements on raw materials, emissions, recycling and power consumption.



Sit-to-stand solutions

Molift's sit-to-stand aids are based on the same ergonomic principles. The products imitate the natural pattern of movement and fulfill the transfer needs of users with some weight bearing capacity.

Molift Raiser Pro

The ergonomic platform for a safe and active sit-to-stand transfer

Early mobilization and rehabilitation

Molift Raiser Pro allows the user to use their own muscle power during a transfer. The natural movement pattern stimulated and the ergonomic design makes the user feel safe during the transfer. Raiser Pro helps the user to rehabilitate in all health environments during every transfer so that mobilization can begin early, already in the intensive care unit.

The stability and maneuverability make Raiser Pro ideal in confined spaces, such as toilets.

Features

- A versatile handle with a soft warm feeling, with an optimal grip for both the user and the caregiver.
- The handle has two different fixing heights for accessories.
- Soft padded legrest with side support, height adjustable with only one hand.
- A central brake mechanism with a one-foot touch that gives audible feedback when activated.
- Low step-on height of 35 mm.

Accessories



Safety strap – Raiser Strap

Provides support around the back or hips when standing up.


S/M Item no. 82508
L/XL Item no. 82509



Safety strap with sliding function – Raiser Strap+

When a user needs to be assisted, Raiser Strap+ can be used by one or two caregivers.

S/M Item no. 82528
L/XL Item no. 82529

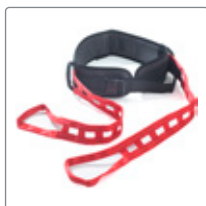

Raiser Belt

Raiser Belt facilitates the work of caregivers who assist users up to a standing position. Can be used by one or two carers

S/M Item no. 82964 (82508 + 82963)
L/XL Item no. 82965 (82509 + 82963)



Product

Molift Raiser Pro

Item no.

16090204

Material

Powder coated aluminum, plastic, steel, chrome 3

Dimensions

Height 1200 mm
Width 540 mm
Depth 630 mm

Weight

Total weight 10,7 kg



For more information visit etac.com





Etac Turner Pro

A sit to stand aid with intuitive design and a small footprint

Safe transfers

Etac Turner Pro is a turn aid with a functional design that offers safe patient turning with standing support. This turn aid enables a good posture and minimises the risk of injury to the carer when transferring. Moreover, it encourages natural standing and sitting movement for the user. The only adjustment necessary is the height of the leg support and the handle. Thereafter, the transfer can be performed safely whilst the caregiver and user can maintain eye contact and communication.


The orange handle provides a strong, high contrasting colour which can help users with impaired vision or dementia. The small footplate makes it the perfect choice of transfer aid when transferring in tight spaces.


Etac Turner Pro is easy to transport due to its low weight of only 7,5 kg, furthermore it is compact and easy to grip.

Features

- Easy to transport due to a low weight of only 7.5 kg. It is compact, evenly balanced and easy to grip.
- The angled leg support is soft and comfortable and provides lateral support. Easy to adjust to whatever height is needed.
- The orange handle provides a strong, high contrasting colour which can help users with impaired vision or dementia.
- The cast aluminium base plate with integrated anti slip function is easy to turn. Works well on different flooring materials including carpets.



Product	Dimensions	Item no.	Material
Etac Turner Pro	800 - 1200 mm	16090106	Cast aluminium, powder coated steel, chrome 3
Dimensions			
Lower leg support height	150 mm		
Lower leg support width	380 mm		
Baseplate diameter	400 mm		
Weight			
Total weight	7.5 kg	For more information visit etac.com	



For more information visit etac.com



Molift Quick Raiser 205

An ergonomic sit to stand hoist with outstanding lifting capacity

Stable, safe and compact

Molift Quick Raiser 205 is a state of the art sit-to-stand hoist, with an electrically adjustable leg base. The hoist is comfortable, stable and safe with exceptional manoeuvrability. The hoisting of Molift Quick Raiser 205 imitates the natural pattern of movement and fulfils the transfer needs of users with some weight bearing capacity.

The unique wheel position, angled push bar and the extraordinary weight distribution of the base, provides a small turning radius. With the small castor option, it fits under most beds and chairs while the electrically adjustable legs enable accessibility almost everywhere. The fixed-leg option benefits a very small footprint.

Features

- Ergonomic design - the inclined column encourages a natural pattern of movement.
- Active 4-point Lifting arm - provides, together with RgoSling Active, extra support when hoisting from sit to stand.
- Easy to operate in moving and handling situations; toilet, sit-to-stand and stand-up training.
- Automatic Warning for periodic inspection and service.



Molift Quick Raiser 205

Model	Battery	Adjustable leg	Fixed leg (FL)
Active 4-point, castor size Ø75mm	SLA	29000	29150
Active 4-point, castor size Ø100mm		29100	29200
StandUp 2-point, castor size Ø75mm		29001	29151
StandUp 2-point castor size Ø100mm		29101	29201
Active 4-point, castor size Ø75mm	Lithium	29002	29152
Active 4-point, castor size Ø100mm		29102	29202
StandUp 2-point, castor size Ø75mm		29003	29153
StandUp 2-point, castor size Ø100mm		29103	29203

For complete variant matrix, contact your local distributor.

Accessories

Accessories	Size	Item no.
StandUp 2-point lifting arm		2910061
StandUp slings		see matrix on next page
Paediatric Footplate		2910065
Set of inlays for footplate		2920151
Leg support, bracket and cushions	Small	2910059
	Large	2910056
Set of cushions	Small	2910076
	Large	2910077
Lithium power pack and charger		2910041

Weight

Total weight (FL)	44 (42,5) kg
-------------------	--------------

Dimensions

External width of wheelbase	620 mm
Internal width of wheelbase	200 mm
Height of footboard:	
Castor size Ø 75/50mm	65 mm
Castor size Ø 100/75mm	90 mm

Technical

Battery SLA	SLA 24V 2,9 Ah
Battery Lithium	LiFePO4 25,6V 2,4 Ah
Charge time	6 hours
Rated performance	100 hoists/75 kg/500 mm
Lifting Range	600 mm
Turning radius	1150 mm
Hoisting speed 75 kg	35 mm/s
Hoisting speed 205 kg	30 mm/s






Molift RgoSling Active

The ideal sling solution for active transfers with Molift Quick Raiser 205


Comfort and easy transfer

Molift Sit to Stand slings are padded and developed for active hoisting together with Molift Quick Raiser 205. The slings provides perfect lumbar support and are the ideal solution for hoisting and moving between different sitting positions. The slings have an efficient anti-slip lining to prevent the sling from sliding up. The waist belt provides additional support and stability.

Molift RgoSling Active is designed for a 4-point lifting arm and has sliding loops for the possibility to an optimal posture during the sit to stand movement. Molift RgoSling Active provides additional grip and support around the hips when standing. The wide back-support distributes weight and pressure evenly over a large surface, which makes it very comfortable for the user.

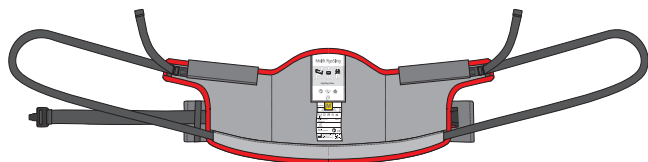

Molift RgoSling Active

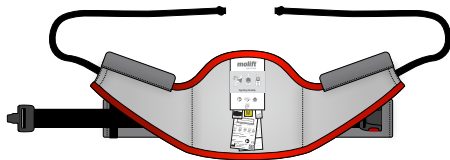

Molift RgoSling StandUp

Molift RgoSling StandUp is designed for a 2-point lifting arm. It is quickly and easily placed around the user. By virtue of the design it is ideal for use in toilet situations.



RgoSling Active and Quick Raiser 205 with 4-point lifting arm

Features

- Perfect solution for hoisting and moving between different sitting positions.
- Self-adjusting sliding loops (RgoSling Active).
- Large and firm support.
- Easy application.

Molift RgoSling

Size

Waist (mm)

Item no.

Active 4-point	XS	500-720	1720910
	S	590-830	1720920
	M	730-950	1720930
	L	960-1210	1720940
	XL	1200-1390	1720950
StandUp 2-point	XS	530-710	1720710
	S	640-850	1720720
	M	740-960	1720730
	L	930-1220	1720740
	XL	1150-1430	1720750
StandUp 2-point with support	XS	530-710	1720810
	S	640-850	1720820
	M	740-960	1720830
	L	930-1220	1720840
	XL	1150-1430	1720850

Information

For more information regarding Active and StandUp lifting arms, see previous..
For complete variant matrix, contact your local distributor.

Material

Polyester



For more information visit etac.com



Mobile Hoists

Molift mobile hoists are well known for many innovative and intuitive solutions such as the 4-point sling bar, integrated service software, a compact lightweight design and an ergonomic, user friendly push handles.



Molift Smart 150

Light, foldable and transportable without tools

Compact smartness

Molift Smart 150 is designed to be easily wheeled and transported to where it is required. It is designed for home care, travelling and easy storage. Molift Smart 150 is foldable in one piece without using any tools.

Due to its low weight, the small footprint and the unique construction, Molift Smart 150 is easy to manoeuvre even in narrow spaces. The hoisting range of 270–1680 mm provides an excellent maximum hoisting height and at the same time gets very low, can lift from the floor. The 4-point sling bar ensures a comfortable and spacious hoisting position for the user.

Molift Smart 150 has been awarded for Design Excellence by the Norwegian Design Council.

Features

- Comfortable and spacious hoisting.
- Compact and foldable without tools.
- Low weight and small footprint.
- Lift from floor.
- Easy to manoeuvre.





Molift Smart 150 folded



Soft travel bag



Travel suitcase

Product		Item no.	Hoisting interval		1410 mm	Material Aluminium, steel
Molift Smart 150 (incl. sling bar)		10200	Width chassis, legs open inner (mm)		925 mm	
			Width chassis, outer (mm)		660 mm	
Accessories		Item no.	Weight		26 kg	 max 150 kg / max 330 lbs 
Soft travel bag		3049177	Total Weight			
Travel suitcase		0990200				
Dimensions			Technical		For more information visit etac.com	
Length		1200 mm	Battery		NiMH 14.4 V – 2.6 Ah	
Leg height		110 mm	Battery charge time		approx. 3 hours	
Dimensions folded		L1160 x W475 x H360 mm	Hoisting speed		60 mm/sec	
Hoisting height		210–1665 mm	Rated performance		40 hoists/75 kg/500 mm	

Molift Partner 255

The flagship of Molift's mobile hoists - with vertical hoisting pattern for excellent stability

Ideal solution for critical care

Molift Partner 255 provides an ideal solution for hoisting in acute and critical care areas in both hospitals and care facilities. The straight vertical column in combination with the fixed sling bar provides excellent stability. It maintains the user at the same distance from the column throughout the hoisting. With an outstanding hoisting range, Molift Partner 255 allows hoisting from the floor as well as high surfaces.

Molift Partner 255 is one of the lightest in its category and combined with excellent construction, ergonomic design and the caster solution, it becomes very easy to manoeuvre even with a lifting capacity of 255 kg.

Molift Partner 255 is available in several variants for the possibility of optimizing the product according to needs; low base chassis, supporting arms for gait/walk training and weight scale.

Features

- High safe working load - 255 kg
- Fixed sling bar and straight vertical column for safe and stable hoisting
- Low weight and easy to manoeuvre
- Suitable for use with Molift stretcher



Product

Molift Partner 255 incl. sling bar

Item no.

13000

For complete variant matrix, contact your local distributor.

Accessories

Casters ø100mm
with directional lock

Item no.

1320140

Weight scale (see page for sling bars)

Sling bars (see page for sling bars)

Stretchers (see page for stretchers)

For more information visit etac.com

Dimensions

Length 1285 mm
Chassis width (outer) 700 mm
Chassis width (inner) 550 mm
Chassis width, open legs (inner) 1065 mm
Hoisting height 390-1680 mm
Hoisting interval 1290 mm
Leg height 115 mm

Weight

Total weight

1285 mm
700 mm
550 mm
1065 mm
390-1680 mm
1290 mm
115 mm

42.5 kg

Technical

Battery
Battery charge time
Hoisting speed
Rated performance

NiMh 26,4 V-2,6 Ah
approx. 3 hours
30 mm/sec
50 hoists/75 kg/500 mm

Material

Aluminium, Steel





Molift Mover

A family of compact, lightweight mobile hoists with great flexibility and large hoisting range. Suitable in homes, hospitals and care environments.

Large hoisting range

Molift Mover has a lightweight aluminium construction that makes the hoist very easy to manoeuvre. The low weight in relation to lifting capacity is the best in class. A large hoisting range allows hoisting from high surfaces as well as hoisting from the floor. The integrated Molift Service Tool ensures correct service and maintenance intervals.

Ergonomic handling

For improved caregiver ergonomics, the Molift Mover hoists are equipped with Molift's unique intuitive push handle. The handle provides multiple and flexible grip possibilities for improved caregiver ergonomics. The easily accessible hand control is very user friendly.

Features

- A lightweight and compact hoist
- Very easy to manoeuvre
- Large hoisting range
- Hoisting from floor

Molift Mover 180

Molift Mover is suitable for transfer situations in home, hospital and care environments. Optional 2- or 4-point sling bar with quick release.



Molift Mover

	Item no.		
Variants (excl. sling bar)	Mover 180	Mover 205	Mover 300
Standard chassis, small casters	18000	14000	15050
Standard chassis, large casters	18003	-	-
Low base chassis	-	14050	-
Standard chassis + supporting arms	-	14150	15150
Medium base chassis + supporting arms	-	-	15151

For complete variant matrix, contact your local distributor.

Accessories

	Item no.
Supporting arms	see etac.com

Accessories

	Item no.
Weight scale	see page for sling bars
Sling bars	see page for sling bars
Stretchers	see page for stretchers
Casters ø100mm with directional lock for Molift Mover 205 and 300	1320140



For more information visit etac.com



Molift Mover 205

Molift Mover 205 is a workhorse for frequent and everyday use, suitable for transfer situations in nursingcare and hospital settings. A quick release function provides optional 2- and 4-point sling bars in different dimensions. A hoist with great flexibility and many accessories.



Molift Mover 300

Molift Mover 300 has the high hoisting capacity of 300 kg. It is the ideal solution for hospitals and care facilities with bariatric clients who have frequent and everyday moving and handling needs. With its low weight and high capacity, it is one of the lightest in its category.



Molift Mover 300 with supporting arm

Facts

for Standard chassis excl. sling bar

	Molift Mover 180	Molift Mover 205	Molift Mover 300
Maximum user weight (SWL)	180 kg	205 kg	300 kg
Total weight	32 kg	37.6 kg	40.1 kg
Length	1260 mm	1285 mm	1285 mm
Chassis width (outer)	700 mm	700 mm	700 mm
Chassis width (inner)	550 mm	550 mm	550 mm
Chassis width, legs open (inner)	1025 mm	1065 mm	1065 mm
Leg height	105 mm	115 mm	115 mm
Battery	NiMH 14.4 - 2.6 Ah	NiMH 26.4 V - 2.6 Ah	NiMH 26.4 V - 2.6 Ah
Battery charge time	3 hrs	3 hrs	3 hrs
Hoisting height	270-1630 mm	2093 mm	2100 mm
Hoisting interval	1365 mm	1431 mm	1420 mm
Hoisting speed	40 mm/sec	40 mm/sec	32.6 mm/sec
Rated performance	50 hoists/75 kg/500 mm	50 hoists/75 kg/500 mm	50 hoists/75 kg/500 mm
Service Software	Included	Included	Included
Material	Aluminium, steel, plastic	Aluminium, steel, plastic	Aluminium, steel, plastic




Sling bars for Molift mobile hoists

Lightweight sling bars with quick release coupling for easy change

Molift aluminium sling bars fit new and older hoists. The lightweight sling bar has a quick-release coupling, which makes it very easy for the caregiver to change.

The weight scale is class III approved and complies to all technical standards. It can calculate BMI (Body mass index) at the touch of a button.

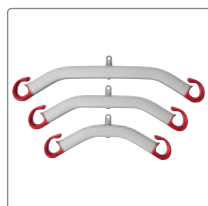

Mover 180 2-point sling bar

Aluminium/zink, SWL 300 kg

Small.....1830043
Medium.....1830042
Large.....1830041



Mover 205 2-point sling bar incl weight scale

Aluminium/zink, SWL 300 kg

Small.....1430120
Medium.....1430121
Large.....1430122



Mover 180 4-point sling bar

Aluminium/zink, SWL 300 kg

Small.....1830110
Medium.....1830111
Large.....1830112
Extra large.....1830113



Mover 205/300 4-point sling bar incl weight scale

Aluminium/zink, SWL 300 kg

Small.....1430120
Medium.....1430121
Large.....1430122



Weight scale Mover 180

SWL 300 kg


Scale.....1840000



Partner 255 4-point sling bar

Aluminium/zink, SWL 255 kg

Medium.....1330160

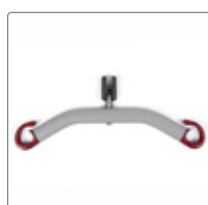

Mover 205 2-point sling bar

Aluminium/zink, SWL 300 kg

Small.....1430117
Medium.....1430116
Large.....1430115



Partner 255 8-point sling bar

Steel, SWL 255 kg

Medium.....2140002



Mover 205/300 4-point sling bar

Aluminium/zink, SWL 300 kg

Medium.....1530101
Large.....1530102
Extra large.....1530103



Weight scale Partner 255

SWL 300 kg

Scale.....1331002

Dimensions (mm)

	Small	Medium	Large	Extra large
2-point Sling bar	H130 x L350	H130 x L450	H130 x L550	-
4-point Sling bar	H210 x W300 x L350	H240 x W300 x L450	H240 x W350 x L550	H240 x W350 x L650
8-point Sling bar	-	H420 x W920 x L1430	-	-

Molift Service Tool


The software that interacts with the hoist

Ensuring correct service and maintenance

All Molift hoists have an integrated service software that indicates when a service is necessary. This ensures user, caregiver and management that the hoist is safe in relation to applicable laws and regulations.

Molift Service Tool calculates the number of lifts and gives a warning signal when a service is required. Through Molift Service Tool, a technician can connect and interact with the hoist to perform and log service and for technical trouble-shooting. The software also allows reading, checking and adjusting key parameters in the hoist as well as logging, reading and printing the service history.

The software provides access to a product specific service checklist.

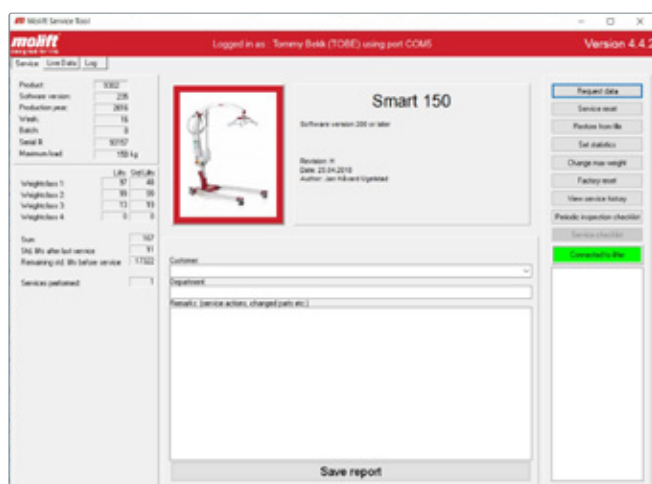

Molift Service Tool enables us to collect data about Molift hoists around the world and to continuously improve our products.

Record and analyse hoist usage

Molift Service Tool can also be used to retrieve data about hoist usage and hoisting patterns. When the data is linked to work-related injury among the employees, it can be a helpful tool for the employers to reduce injuries. Furthermore Molift Service Tool can improve efficiency as it indicates whether the types of hoists in use are suitable for the patient group and their needs.





Ceiling Hoists

Molift Ceiling Hoists and Molift Rail System are the optimal solutions to assist caregivers. The hoists are intuitive and safe to use and can solve all hoisting needs for disabled people and improve accessibility to the user and increases caregiver and user safety.

Molift Nomad

The portable low weight hoist, suitable for all day-to-day operations.

Low weight - high hoisting capacity

The weight of the Molift Nomad is only 6.8 kg. Despite the low weight it has a hoisting capacity up to 255 kg. Molift Nomad has a unique curved design and the integrated 4-point sling bar, used with any Molift sling, provides the user with a comfortable and spacious position. Molift Nomad is equipped with numerous safety features such as electronic overload protection, electrical emergency lowering and emergency stop.

Easy room-to-room transfer

Molift Nomad has integrated carry handles for easy handling. The connection is easy and no tools are required. In combination with the special designed trolley, Molift Nomad can be used to climb from one rail system to another without a transition rail. This results in easy transfer from room to room. By using Molift Nomad Extension Arm the hoist can easily be connected and disconnected to the trolley from the floor. Molift Nomad Cart is ideal for storing and transporting the hoist.



Features

- Low weight and high hoisting capacity.
- Portable ceiling hoist.
- Climbing room to room.
- Easy connection without tools to trolley.
- Integrated 4-point sling bar.

Description	Item no.	Dimensions	Weight	Technical	Hoisting speed	Rated performance	Hoisting range	Material
Molift Nomad 205 (SWL 205 kg)	22205	Length	522 mm		30 mm/sec			
Molift Nomad 255 (SWL 255 kg)	22255	Width	200 mm		80 hoists/75 kg/500 mm			
		Height	330 mm		Up to 2000 mm			
Accessories	Item no.	Weight	Technical					
Weight scale, class III	1100253	Total weight	6.8 kg					
Cart	1200005							
Extension arm L:900 mm	1200004							
Extension arm L:1050 mm	1200006							
				NiMh 26.4 V – 2.6 Ah approx. 3 hours				



For more information visit etac.com





Molift Air

For a better working situation for carers and a safe and dignified quality of care for users

Dignified quality of care

Molift Air is a smart and user-friendly family of ceiling hoists, with outstanding lifting capacity. Molift Air is suitable for sitting and horizontal transfers as well as standing and gait training situations.

Together with Molift Rail System and a wide range of accessories, Air ceiling hoists offers great flexibility and eliminates strenuous manual lifting in a variety of care situations. A complete solution within patient handling, that creates a better working situation for carers and provides a better and more dignified quality of care for users.

Light and flexible

Molift Air is used with fixed ceiling rails or portable free-standing systems. The hoist is easy to remove and replace during installation and servicing, due to very low weight and the clever quick-release function. Together with all safety features, a large hoisting range and the multiple sling bar range, is this the complete solution in all health care settings.

Molift Air can easily be adapted to other manufacturers rail systems on the market.



Features

- Small hoist - great lifting capacity
- Flexible solutions
- Low lifecycle cost
- Silent and soft operation.

Molift Air 200

Molift Air 200, with its lifting capacity of 200 kg, is optimized for long term care facilities, nursing homes and private properties. The hoist is charged through the handset.

Product description	Battery	Item no.
Molift Air 200	SLA	26200
Molift Air 200	Litium	22255
Molift Air 350	SLA	26350
Molift Air 350	Litium	26351
Molift Air 350 Propulsion	SLA	26352
Molift Air 350 IRC and Propulsion	SLA	26353
Molift Air 350 IRC	SLA	26354
Molift Air 500	SLA	26500

Accessories	Item no.
Sling bars	see page for sling bars
Weight scale	see page for sling bars
Room to room transfer unit	2510295

Material

Aluminium, plastic



For more information visit etac.com



Molift Air 350

Molift Air 350 has an outstanding lifting capacity and is optimized for long lifetime and low service costs. Molift Air 350 has several variants available; In-Rail charging, propulsion trolley and lithium battery cells.

Molift Air 500

Molift Air 500 is designed for bariatric users and has a lifting capacity of 500 kg. With the specially designed sling bar, the bariatric user can be lifted safely and securely.



Facts

	Molift Air 200	Molift Air 350	Molift Air 500
Maximum user weight (SWL)	205 kg	350 kg	500 kg
Total weight (with SLA battery)	8.6 kg	8.6 kg	15 kg
Length	360 mm	360 mm	360 mm
Width	190 mm	190 mm	250 mm
Height	195 mm	195 mm	195 mm
Battery	26.4 V SLA 2.6 Ah	26.4 V SLA 2.6 Ah	26.4 V SLA 2.6 Ah
Battery charger	100–240 V AC, 40–60 Hz	100–240 V AC, 40–60 Hz	100–240 V AC, 40–60 Hz
Battery charge time	6 hrs	6 hrs	6 hrs
Motor	24 V DC	24 V DC	24 V DC
Hoisting interval	3000 mm	2600 mm	3000 mm
Hoisting speed 75 kg	60 mm/sec	60 mm/sec	30 mm/sec
Expected lifetime	30 000 lifts SWL/10 years	30 000 liftsw SWL/10 years	30 000 lifts SWL/10 years
Molift Service Software	Included	Included	Included
Emergency lowering	Mechanical and electrical	Mechanical and electrical	Mechanical and electrical
Emergency stop	Reachable from floor with string. Reset possible from floor		
Protection class; hand control & hoist motor	IP24	IP24	IP24



Molift Rail System

A cost-effective and flexible solution for all environments

Efficient installation process

Molift Rail System has rail profiles with innovative solutions for integration into all types of ceilings and walls. The wide range of fixing options results in making it one of the most most efficient installation processes in the market. It is easily integrated and is aesthetically adjusted to the design and architecture of the room. The range of switches and other transition options makes the rail system extremely flexible.

Customised capacity

Depending on the different needs and choice of the hoist motor, the Molift Rail System can be constructed from 160 up to 500 kg.

Cost effective solutions

The efficient installation process, the flexibility and the many features make Molift Rail System a cost-effective, competitive and profitable investment for hospitals, institutions, health care settings and the home. This helps create a practical and efficient workplace for nurses and caregivers, as well as an attractive environment for patients and visitors.

Traverse rail system

One of the most flexible solutions on the market, covering all needs of hoisting, regardless of renewal or equipment changes in the room.

Traverse switch

Enables two different rail systems to be connected, for example between a traverse system and a single rail. Very useful for transfers from one room to another.

Single rail system


Straight rail and/or curve rails gives a simple and effective patient transfer for the most common lifting situations such as bed-to-wheelchair hoisting.

Room to room transfer

When transfer is needed from one room to another or from one rail system to the next system, Molift can provide Room to Room function (Climbing), switch systems or continuing rail systems through doorways etc.

Rail switch

Rail Switch - enables a change of direction of single rails, for example from a straight rail to a curve rail.

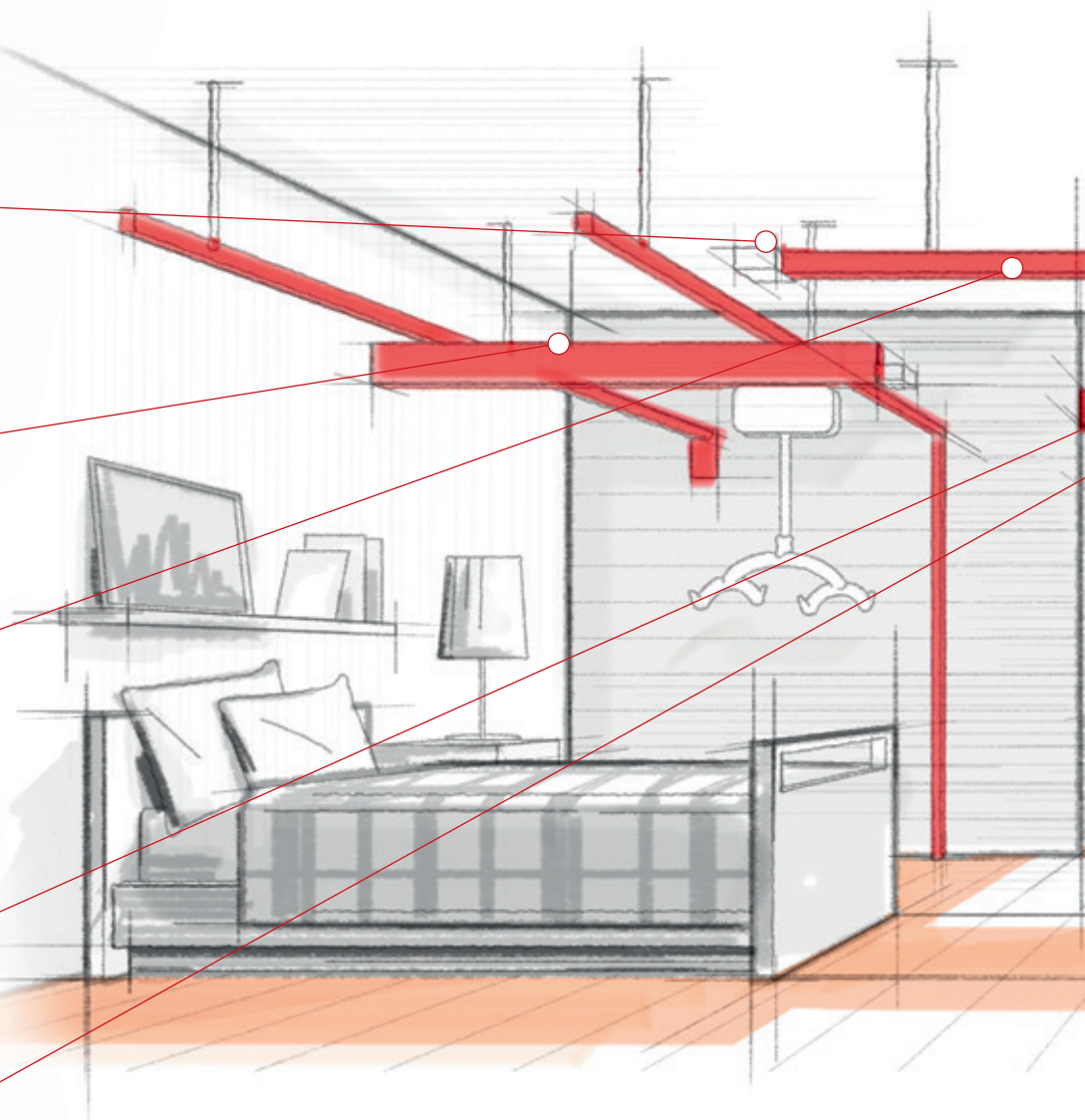

Flexible system

- Easy to install.
- Various weight classes.
- A variety of installation options – integrated in, or under the false ceiling, on the wall, on the floor with upright supports, etc.
- Can be installed in all types of ceilings, walls and floors.

Possibilities

Straight profiles
Curves - 30, 45, 60 and 90 degrees
Fixed brackets, telescope brackets
Flexible traverse options
Traverse switch
Rail switch
In-rail charging
Curtain pass through

Maximum user weight (SWL)

160–500 kg

Colours

White

Rails


3 dimensions:
Height: 62, 112, 142 mm

Max distance between fixtures

Up to 6000 mm

Switches

Rail switch
Traverse switch

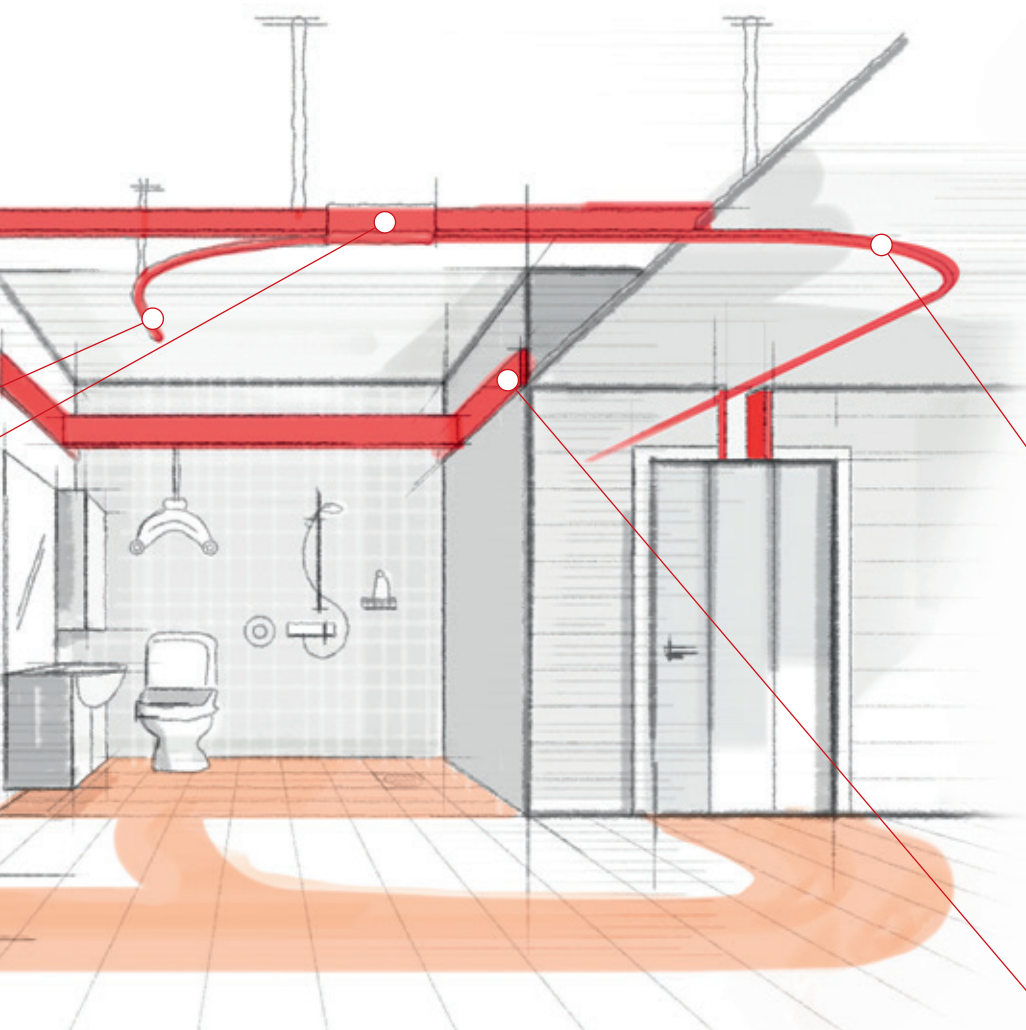

In-rail charging

In rail charging allows the hoist to charge batteries automatically anywhere in the system. This way the hoist is always ready for operation. In-rail charging can be mounted in straight rail and traverse rail systems only.

Integrated rail

Integrated Rail System - a single integrated rail offers a very discrete and clean solution.

Curtain pass through

Curtain pass through allows gaps for curtains mounted in the ceiling of the ward. This means several beds can share one motor, and still maintain privacy when individual nursing is required.

Wall mounted rails

Wall mounted rails - if the ceiling cannot support the system, it is possible to mount the rail on the wall.



Sling bars for Molift ceiling hoists

Lightweight sling bars with quick release coupling for easy change

Molift aluminium sling bars fit new and older hoists. The lightweight sling bar has a quick-release coupling, which makes it very easy for the caregiver to change.

The weight scale is class III approved and complies to all technical standards. It can calculate BMI (Body mass index) at the touch of a button.



Molift Air 2-point sling bar

Aluminium/zink, SWL 350 kg

Small.....2530115
Medium.....2530120
Large.....2530125



Molift Air 4-point bariatric pigtail sling bar

Steel, SWL 550 kg

Large.....2510435



Molift Air 2-point pigtail sling bar

Aluminium, SWL 350 kg

Small.....1830033
Medium.....1830032
Large.....1830031



Molift Air 8-point sling bar

Steel, SWL 300 kg

Medium.....2140003



Molift Air 4-point sling bar

Aluminium/zink, SWL 350 kg


Small.....2530210
Medium.....2530211
Large.....2530212
Extra large.....2530213



Weight scale Molift Air

(not available for 8-pt sling bar)

Scale SWL 300 kg.....1840000
Scale SWL 400 kg.....1431021

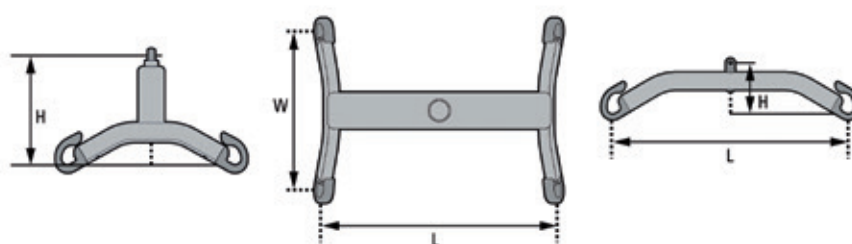

Dimensions (mm)

	Small	Medium	Large	Extra large
2-point Sling bar	H130 x L350	H130 x L450	H130 x L550	-
2-point Sling bar pig tail	H150 x L400	H150 x L500	H150 x L600	-
4-point Sling bar	H210 x W300 x L350	H240 x W300 x L450	H240 x W350 x L550	H240 x W350 x L650
4-point Sling bar pig tail	-	-	H200 x W550 x L 550	-
8-point Sling bar	-	H420 x W920 x L1430	-	-

Molift Quattro

The temporary hoist solution with excellent stability

A free-standing traverse gantry

Molift Quattro is a free-standing hoist unit that offers a (four post) Traverse System which is the best solution for temporary hoisting needs or when the facility is not built for fixed installations. The gantry is easy to install without the need to mount additional fixtures on walls or ceilings. It is ideally suited for temporary hoist applications in nursing or private homes.

For overhead lifting situations

The traverse has a lightweight aluminium frame and four steel posts that provides excellent stability and covers the surface of the hoisting requirement. Molift Quattro solves all overhead lifting situations together with the Molift Air or Molift Nomad hoist motor and a sling. With a lifting capacity of max 300 kg, it is ideal for lifting and transferring a person from or to bed, floor, chair, wheelchair, toilet etc.

Features

- Easy to install - designed for quick installation by two persons (preferably).
- No wall or ceiling mounting is required.
- Supreme stability and manoeuvrability .
- Easy to charge - the hoist motor can easily be charged through the hand control or through the in-rail charging option.



Description

Molift Quattro 3	1600030
Molift Quattro 4	1600035

Accessories

Adjustable foot (set of 4)	1610200
----------------------------	---------

Dimensions

Boom lengths	3000/4000 mm
Height	2330 mm
Height with adjustable foots	2360 - 2430 mm

Internal height	2190 mm
Width/Depth	3160 - 4160 mm
Maximum internal space	3000/4000 x 2190 mm
The length of the posts are fixed. The lowest allowed distance from floor to ceiling is 2400 mm.	
Molift ceiling hoists, Air and Nomad are ordered separately	

Weight

Total weight Quattro 3 (complete unit)	111 kg
Total weight Quattro 4 (complete unit)	132 kg

Technical

Expected Lifetime:
30 000 hoisting cycles (SWL) or 10 years

Material

Steel, aluminium, plastic



For more information visit etac.com



Molift Duo

A flexible and free standing gantry system for temporary hoist situations

The easy solution for most room sizes

Molift Duo offers a system that is easy to install without being fixed to the ceiling or walls. It is ideally suited to temporary hoist applications in hospitals, nursing homes or private homes. In combination with Molift Air or Molift Nomad hoist and a sling, Molift Duo solves most overhead lifting situations, whatever the environment. With a lifting capacity of 300 kg, it is ideal for lifting and transferring a person from or to; bed, floor, chair, wheelchair, toilet etc.

Lightweight construction

Molift Duo has a lightweight aluminium construction and a functional design which makes the gantry system easy to install, dismantle and move elsewhere.

This 2-post gantry system has a rail length between 2000 mm and 3500 mm and has 4 step adjustable heights up to 2550 mm, making this free standing hoist unit suitable for all bed and room sizes.

Features

- Easy to install - due to the lightweight system, it is easy to install the gantry and move it to a new location.
- Supreme stability - the gantry system has supreme stability and a functional design.
- Adjustable height - height positions can be set at four fixed positions from 225 cm up to 255 cm. Adjust the height with the quick release pin.
- Charging of lift motor - Molift hoist motors can easily be charged thorough the hand control or optionally through in-rail charging.



Description

Molift Duo 2000	1600013
Molift Duo 2500	1600015
Molift Duo 3000	1600017
Molift Duo 3500	1600019

Item no.

Dimensions

Rail lengths	2000/2500/3000/3500 mm
Height	2250 - 2550 mm
Depth	1200 mm
Maximum internal space	3500 x 2400 mm

Weight

Total weight (complete unit)	38 - 44,3 kg
------------------------------	--------------

Material

Anodized aluminium, Plastic



For more information visit etac.com





Stretchers

Molift offers a range of accessories for different needs. One of them is the unique Molift Strecher which can be used directly on a standard Molift 4-point sling bar.



Molift RgoSling Fabric Stretcher

A soft fabric hoisting device for horizontal transfers at hospitals and institutions

The soft stretcher option

Molift RgoSling Fabric Stretcher is a soft fabric hoisting device for hospitals and institutions. It is designed for horizontal transfers. Molift RgoSling Fabric Stretcher is suitable when hoisting in a normal sling would be inappropriate due to medical conditions, for example comatose or sedated clients that require a stable recumbent position. It is a comfortable alternative to a stretcher. Molift RgoSling Fabric Stretcher is compatible with Molift Partner 255 and Molift Air.

8-point sling bar

Hoisting with this product requires an 8-point sling bar. The sling bar is compatible with Molift Partner 255 and Molift Air. Molift Fabric Stretcher has integrated hoisting straps that are hooked directly to the sling bar, with some possibilities to adjust in length by using the loops. The multiple hoisting points provide an open and stable horizontal transfer.

Features

- Molift RgoSling Fabric Stretcher is suitable when hoisting in a normal sling would be inappropriate due to medical conditions.
- The application is easily performed, with an Immedia SlingOn, when the client is in a horizontal position.
- Open and stable horizontal hoisting - the 8-point sling bar with multiple hoisting points provides an open and stable horizontal transfer.



Description

Description	Item no.
Molift RgoSling Fabric Stretcher M	1721520
Molift RgoSling Fabric Stretcher L	1721530
Molift RgoSling Fabric Stretcher XL	1721540

Dimensions

Length outer	Medium	1750 mm
	Large	1780 mm
	Extra large	1830 mm
Width outer	Medium	450 mm
	Large	650 mm
	Extra large	900 mm

Material

Molift Fabric Stretch:

Polyester



* The safe working load will always dependent on the hoists safe working load when less than the Fabric Stretcher.

For more information visit etac.com

Molift Stretcher

A scoop stretcher, ideal for horizontal hoisting situations

Easy to use

Molift Stretcher is a scoop stretcher for hospitals and institutions. It is designed for horizontal transfers and when sling hoisting is considered inappropriate due to medical conditions. Molift Stretcher is an optimal solution for transfer of clients who are comatose or sedated, as well as clients with suspected spinal injuries and/or multiple fractures.

Ideal for emergency care

In combination with Molift standard 4-point sling bar, it is ideal for emergency care needs, when time is critical and the patient has limited mobility. Molift Stretcher can be divided in two parts in length and is therefore easy to place under the client.

Molift Stretcher has five handles along each side, which provides a manual hoisting alternative from the floor and allows horizontal and lateral transfer between beds and surfaces.

Features

- A unique patented scoop design, for safe and efficient moving.
- Low total weight and easy to carry.
- A safe locking system - easily opened when unloaded - completely locked when loaded.
- Many possible uses and compatible with X-ray, CT and MRI. Designed for minimal shadow in X-ray.



Description

Molift Stretcher MRI

Item no.

2150108

Dimensions

Height
Length outer
Width outer

126 mm
1980 mm
630 mm

60-85° max 300 kg max 661 lbs



Material

Plastic composite body
Stainless steel locks

Compatibility

Molift Stretcher MRI:
Compatibility with MRI. Approved for X-ray and CT-scan

* The SWL always dependent on the hoists SWL when less than the stretcher.

For more information visit etac.com

Weight

Total weight

14.5 kg



Molift RgoSling 4-point

Molift RgoSlings offer a wide range of models and sizes ensuring we have a solution, whatever the users requirements. Molift RgoSlings are designed for a 4-point sling bar.

Sliding loops

Sliding loops on head straps. Making a smooth transfer from lying to sitting position. Ensures a comfortable and safe position for the patient.

Colour coding

Colour coding of hoisting loops helps you find the right fit.

Comfortable padding

Padding in back and leg support distributes the weight over a larger area, making the sling more comfortable for the user.

Labelling

Clear, practical information on a highly visible label. Washing instructions and technical information are found in a protective pocket at the back.

Guide straps

Guide straps assist with positioning during the lifting process.

Application pocket

Pockets in the back and on the leg supports make application easier.



Molift RgoSlings

Molift has a long reputation for manufacturing high quality hoists and slings. RgoSling is a portfolio of 4-point slings developed to meet the demands for the everyday needs in healthcare.

Standard models: HighBack and MediumBack

Our standard models are designed for use in the usual transfer situations and are suitable for most users. The standard models are available in two back heights: HighBack and MediumBack.

Models with specific features

We know that users have many different requirements, so we've designed additional models with specific features to meet their needs, including the RgoSling Toilet and RgoSling Amputee.



Molift RgoSling HighBack

Standard leg support opening

- Sliding loops on head straps.
- Support for the entire back and the head.
- Hoisting to and from a lying or a sitting position.
- All-round sling that fits most users.
- Hoisting from floor.

Molift RgoSling MediumBack

Standard leg support opening

- Sliding loops on head straps.
- Suitable for a user with headcontrol.
- All-round sling that fits most users.
- Adjust the sitting angle by the leg straps.

Size		HighBack padded	HighBack net	Size		MediumBack padded	MediumBack net
XXS		1720200	1720400	XXS		1720100	1720300
XS		1720210	1720410	XS		1720110	1720310
S		1720220	1720420	S		1720120	1720320
M		1720230	1720430	M		1720130	1720330
L		1720240	1720440	L		1720140	1720340
XL		1720250	1720450	XL		1720150	1720350
XXL		1720260	1720460	XXL		1720160	1720360

Maximum user mass: 300 kg
HighBack padded XL - XXL 350 kg

Maximum user mass: 300 kg
MediumBack padded XL - XXL 350 kg



Molift RgoSling Toilet

Wide leg support opening

- Fixed loops.
- Easy access sling, facilitate dressing and undressing.
- For users with good core stability and muscle tonus.
- Comfortable padding.
- Flexible leg support.
- Belt included - for extra security.

Molift RgoSling Amputee

Narrow leg support opening

- Sliding loops on head straps.
- For users with an amputated leg or a user with reduced core stability or muscle tonus.
- Ideal solution for user with a high amputee.
- Comfortable padding.
- Belt included - for extra security.

Size		LowBack	HighBack
XXS		1720500	1720600
XS		1720510	1720610
S		1720520	1720620
M		1720530	1720630
L		1720540	1720640
XL		1720550	1721450
XXL		1720560	1720660

Size		MediumBack	HighBack
XXS		1721100	1721200
XS		1721110	1721210
S		1721120	1721220
M		1721130	1721230
L		1721140	1721240
XL		1721150	1721250
XXL		1721160	1721260



Molift RgoSling Shadow

Standard leg support opening

- Sliding loops on head straps.
- Suitable for users in moulded seats.
- No straps, seams, pockets or labels in contact with the user.
- Polyester net with thin webbing: no straps, seams, pockets or labels in the contact surface to the user.

Molift RgoSling Comfort HighBack

Full body support sling without split leg

- Full body and head support.
- Can be left in the wheelchair underneath the user.
- Polyester net which provides good ventilation.
- Red seam indicating the center line of the back, to facilitate correct application.

Size		HighBack	Size		HighBack
XXS		1722900	XS		1721410
XS		1722910	S		1721420
S		1722920	M		1721430
M		1722930	L		1721440
L		1722940	XL		1721450
XL		1722950	XXL		1721460



Molift RgoSling Ambulating Vest

For gait training with assistance

- 4 straps with fixed loops.
- Self tightening strap in the waist.
- Standing or balance exercise.
- Double buckles for extra security.
- Standard polyester fabric with padding.
- Padded groin straps is ordered separately.

Molift RgoSling Fabric Stretcher

A soft fabric hoisting device for horizontal transfers



- Molift RgoSling Fabric Stretcher is suitable when hoisting in a normal sling would be inappropriate due to medical conditions.
- Application easily performed.
- Standard polyester fabric with padding of the head part.
- Used with an 8-point sling bar.

Size		Item no
XXS		1721300
XS		1721310
S		1721320
M		1721330
L		1721340
XL		1721350

Size	Groin straps	Item no
XXS-XS	single pack (order two for a pair)	1721700
S-M		1721720
L-XL		1721740

Size	Width x Length	Item no
M	450 x 1750 mm	1721520
L	650 x 1780 mm	1721530
XL	900 x 1830 mm	1721540



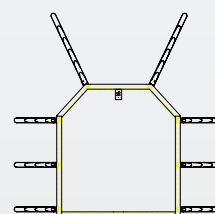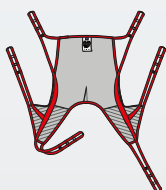

Molift RgoSling MediumBack Plus

Full body support for users weighing from 160 kg

- Provides full body support
- Straps with fixed loops and multiple height options
- Comfortable padded leg-support with dual straps
- Colour coded straps
- Split leg opening
- Can be used with a separate headrest

Molift RgoSling Repositioning sheet

An in-bed lift and positioning sheet

- Turning and transferring patients higher in the bed.
- Lifting from bed to another supine surface, stretcher etc.
- Stays in bed - for comfortable and dignified care.
- Facilitates care situations - saves time for caregivers.
- Colour coded straps facilitate communication when attaching the straps equally on both sides.

Size

L
XL



MediumBack Plus

1721940
1721950

Maximum user mass: 500 kg

Size

L 1340 x W 1300 mm
L 1340 x W 1300 mm

Quantity

1 pc
10 pcs

Repositioning sheet

1721570
1721570-C

Maximum user mass: 500 kg



Molift EvoSling 2-point

Molift RgoSling offer a wide range of sling models and sizes ensuring we have a solution, whatever the users requirements. Molift EvoSlings are designed for a 2-point sling bar.

Adjust with fixed loops

Hoisting loops on both head and leg supports allow the distance to the sling bar and the user's sitting angle to be adjusted.

Colour coding

Colour coding of hoisting loops helps you find the right fit.

Comfortable padding

Padding in back and leg support distributes the weight over a larger area, giving a much more comfortable lift for the user.

Labelling

Clear, practical information on a highly visible label. Washing instructions and technical information are found in a protective pocket at the back.

Guide strap

Guide straps assist with positioning during the lifting process.

Application pocket

Pockets in the back and on the leg supports make application easier.



Molift EvoSlings


Molift has a long reputation for manufacturing high quality hoists and slings. EvoSling is a portfolio of 2-point slings developed to meet the demands for the everyday needs in healthcare.

Standard models: HighBack, MediumBack & LowBack


Our standard models are designed for use in the usual transfer situations and are suitable for most users. They are available in three different back heights: HighBack, MediumBack and LowBack. The padded version has polyester fabric with padded back and legs and net version has polyester net with padded legrests.

Models with specific features


We know that users have many different requirements, so we've designed additional models with specific features to meet their needs, including the EvoSling Hygiene and EvoSling Ampu.

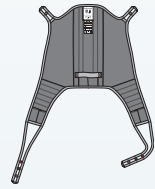


HighBack padded

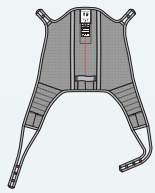

HighBack net

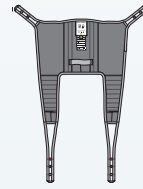

MediumBack padded

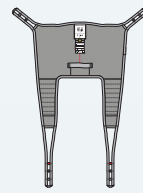

MediumBack net



Molift EvoSling HighBack

Standard leg support opening

- Support for the entire back and the head.
- Hoisting to and from lying or sitting position.
- Hoisting from the floor.
- All-round sling that fits most users.
- Net version has padded legs for extra comfort.

Molift EvoSling MediumBack

Standard leg support opening

- Suitable for a user with headcontrol.
- All-round sling that fits most users.
- Net version has padded legs for extra comfort.

Size		HighBack padded	HighBack net padded	Size		MediumBack padded	MediumBack net padded
XS	<input type="checkbox"/>	1722810	1722710	XS	<input type="checkbox"/>	1723010	1723110
S	<input checked="" type="checkbox"/>	1722820	1722720	S	<input checked="" type="checkbox"/>	1723020	1723120
M	<input type="checkbox"/>	1722830	1722730	M	<input type="checkbox"/>	1723030	1723130
L	<input checked="" type="checkbox"/>	1722840	1722740	L	<input checked="" type="checkbox"/>	1723040	1723140
XL	<input type="checkbox"/>	1722850	1722750	XL	<input type="checkbox"/>	1723050	1723150
XXL	<input type="checkbox"/>	1722860	1722760	XXL	<input type="checkbox"/>	1723060	1723160



Molift EvoSling LowBack

Standard leg support opening

- Suitable for users with head control.
- Comfortable reclined sitting position.
- All-round sling that fits most users.
- Net version has padded legs for extra comfort.


Molift EvoSling Hygiene

Wide leg support opening

- Comfortable padding.
- Belt for extra security.
- For users with good core stability and muscle tonus.
- Easy access sling, facilitate dressing and undressing.

Size		LowBack padded	LowBack net padded
XS	<input type="checkbox"/>	1722310	1722510
S	<input checked="" type="checkbox"/>	1722320	1722520
M	<input checked="" type="checkbox"/>	1722330	1722530
L	<input checked="" type="checkbox"/>	1722340	1722540
XL	<input checked="" type="checkbox"/>	1722350	1722550
XXL	<input type="checkbox"/>	1722360	1722560

Size		LowBack padded
XS	<input type="checkbox"/>	1722210
S	<input checked="" type="checkbox"/>	1722220
M	<input checked="" type="checkbox"/>	1722230
L	<input checked="" type="checkbox"/>	1722240
XL	<input checked="" type="checkbox"/>	1722250
XXL	<input type="checkbox"/>	1722260

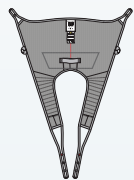

MediumBack
net padded



HighBack net



Molift EvoSling Ampu MediumBack

Narrow leg support opening

- For users with an amputated leg or a user with reduced core stability or muscle tonus.
- For user with headcontrol.
- Anti-slip material on part of the inside.
- Reclined sitting position.
- Comfortable leg padding.
- Polyester net.

Molift EvoSling Shadow HighBack

Standard leg support opening

- For users who need extra support.
- Suitable for user in moulded seats.
- Six straps gives more support, security and flexibility.
- Arms can be placed inside or outside the sling.
- Polyester net with thin webbing: no straps, seams, pockets or labels in the contact surface to the user.

Size		MediumBack net padded	Size	HighBack
XS		1723210	XXS	1723900
S		1723220	XS	1723910
M		1723230	S	1723920
L		1723240	M	1723930
XL		1723250	L	1723940
XXL		1723260	XL	1723950



Molift EvoSling Comfort MediumBack

Full body support sling without split leg

- Ideal for users with joint problems or sensitivity to pain.
- Can be left in the wheelchair underneath the user.
- Polyester net, which provides good ventilation.
- Red seam indicating the center line of the back, to facilitate correct application.

Molift EvoSling FlexiStrap

Avoids heavy lifting of a leg or an arm

- Lifts a leg or arm to avoid heavy lifting.
- The FlexiStrap available in two widths.
- Padded polyester for comfortable lifting.
- Pocket on one side to facilitate application.

Size

XS	<input type="checkbox"/>
S	<input checked="" type="checkbox"/>
M	<input type="checkbox"/>
L	<input type="checkbox"/>
XL	<input type="checkbox"/>
XXL	<input type="checkbox"/>

MediumBack net

1722610
1722620
1722630
1722640
1722650
1722660

Version

Molift EvoSling FlexiStrap, standard
Molift EvoSling FlexiStrap, wide

Width x Length

15 x 95 cm
25 x 95 cm

FlexiStrap

1723310
1723330



Molift UnoSling

Slings for short-term use are specific to a single user, a patient specific sling. An effective solution for preventing the transmission of infection.

Molift UnoSling are designed to fit a 2-point or a 4-point sling bar.



Molift UnoSling HighBack

Standard leg support opening

- A personal disposable sling for short-term use.
- Hoisting to and from lying or sitting position.
- Support for the entire back and the head.
- For environments with a short term need for a sling,. Can be used several times by the same user.



Disposable product supplied in packaging which makes it easy to store and take out a new sling when needed.

Size	Number in package	Art. no.
S	20	3060020
M	20	3060030
L	20	3060040
XL	20	3060050
XXL	10	3060060

Disposable product

Do not wash. The UnoSling should be discarded when the user no longer needs it, or if it becomes damaged or dirty.



max
300
kg


max
660
lbs






Facts

Materials

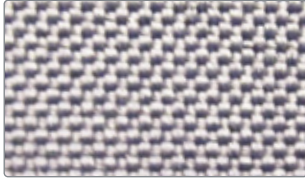

Polyester

The standard polyester fabric is durable and hard wearing. The material is easy to maintain and has a dirt repellent surface.

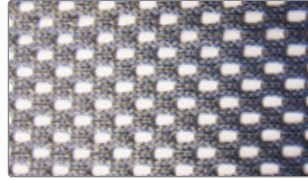

Polyester net

Hard wearing and durable soft net fabric, breathes well and dries quickly. Suitable for shower and bathing situations.



Non-woven

Polypropylene non-woven. Patient specific slings when the sling is for short term use.



Durable polyester

Extra strong durable polyester fabric for products that has a weight limit up to 500 kg. The material is easy to maintain and dries quickly.

Detailed information

Size	Recommended user weight
XXS	12 - 17 kg
XS	17 - 25 kg
S	25 - 50 kg
M	45 - 95 kg
L	90 - 160 kg
XL	160 - 240 kg
XXL	230 - 300 kg

Washing instructions



Washing temperature 60-85 °C / 140-185 °F
Washing at higher temperatures wears out the material faster.



All slings have a maximum user mass of 300 kg.
Some exceptions exist which are specified under each product.



Molift Slings accessories




Molift EvoSling Support

Accessories for slings without head support
The head support is attached to the sling and its length can be adjusted as required using the plastic buckle. It gives the user support for the head and is a useful addition when the user requires more support.



Molift EvoSling Soft Padding

For pressure-sensitive users
EvoSling Soft Padding consists of pressure-relieving padded sleeves which are slipped over the leg parts of slings with split-leg support. Washing 40 °C / 104 °F.

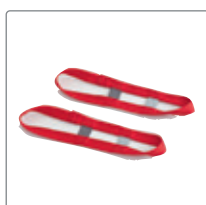

Molift RgoSling Extension Loops

Straps for easier adjustment
RgoSling: Used to extend the hoisting straps, flexible for 3 different length.
N.B: Order 2 for a pair.



Molift EvoSling Extension Loops

Straps for easier adjustment
Used to extend the hoisting straps or to keep the leg supports together.
N.B: Order 2 for a pair.



Molift Strap Marker

For marking a selected hoisting loop for head or leg support hoisting straps. Mark the hoisting loops selected during individual testing of the sling.



Immedia SlingOn

When the user can't be turned onto their side
Immedia SlingOn used with a folding technique makes it possible to apply the sling without first turning the user over to the side.
Washing in 80 °C or 176 °F.

Product		item no.
Molift EvoSling Support		1723400
Product (sold in pairs)	Size	item no.
Molift EvoSling Soft Padding	XS - M	1723710
Molift EvoSling Soft Padding	L - XXL	1723740
Molift EvoSling Hygiene Soft Padding	S - XXL	1723810
Product	Extension	item no.
Molift RgoSling Extension loops (flexible)	+60, 180 or 260 mm	1721600

Product	Extension	item no.
Molift EvoSling Extension loop	+90 mm	1723510
Molift EvoSling Extension loop	+180 mm	1723520
Product (sold in pairs)		item no.
Molift Strap Marker		1723600
Product (sold in pairs)	Size	item no.
Immedia SlingOn - Small	800 x 1000 mm	IM437/S
Immedia SlingOn - Medium	800 x 1700 mm	IM438
Immedia SlingOn - Large	1200 x 1500 mm	IM437



Etac is a world-leading developer and provider of ergonomic assistive devices and patient handling equipment. Our heart lies in the solutions that optimise quality of life for the individual, their family and caregivers.

For the latest news and continuously updated product information – visit: www.etac.com