

Mobility Designed Forearm Comfort Crutch

HCPCS:



FIGURE A



FIGURE B



FIGURE C



FIGURE D

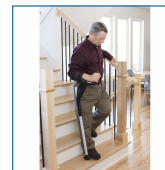


FIGURE E



FIGURE F



FIGURE G

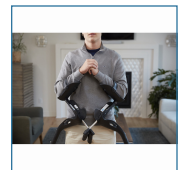


FIGURE H

- A hinged arm cradle that can be unlocked to provide an enhanced range of motion enabling the user to reach for things without having to remove the crutch.
- Adjustable, ergonomic hand grips rotate out of the way with the push of a button, freeing your hands for other things.
- Flexible armbands made of anti-microbial material that conform to each users arm, providing support while allowing the user to easily get in and out of the crutches.
- Removable, shock absorbing feet help users to navigate a variety of surfaces. The patented design allows water to disperse through hollow channels much like a car tire, to avoid slippage on wet surfaces.
- Fully adjustable in both arm length and height accommodating users 411 to 68
- Handle Length: 11"to 17"
- MD20041: Box includes Two Fully Assembled (Tool-Free) Pairs of Crutches with complete set of Training Support to help maximize usage.
- MD30002: Box includes Three Non-Assembled Pairs of Crutches with complete set of Training Support to help maximize usage.
- Blue button - Right Side
- Black button - Left Side
- The patent pending Comfort Crutch's unique design cradles your elbows to evenly distribute your weight, allowing the pressure to be removed from your wrists, armpits and hands.
- The M+D Crutch offers a wide variety of one of a kind feature including: a hinged arm cradle that can be unlocked to provide enhanced range of motion, customizable braces with flexible, quick release straps that keep arms in place when they need to be and free when they don't, shock absorbing feet to help users navigate a variety of surfaces and rotating handles allow for more mobility!
- MAXIMUM MOVEMENT:
- HANDS-FREE:
- QUICK RELEASE STRAPS:
- SHOCK ABSORBING FEET:
- CUSTOMIZABLE FIT:
- SIMPLE SETUP:
- Ergonomically Designed to fit either your left or right hand, the M+D Crutch is color coded and has a braille dot so you can easily know which crutch is which!
- Aluminum tubing features an oval shape providing a lighter weight product front to back

PRODUCTS	Description
MD20041	Fully Assembled MD Crutch, Blk 2PR/BX
MD30002	Non-Assembled MD Crutch, Blk 1PR,3PR/BX

SPECIFICATIONS	Description
Contains Latex	No

SPECIFICATIONS	Description
Height Adjustable	39" - 52"
Material	Aluminum
Product Weight	2.9 lbs per crutch
Retail Box Dimensions	13" (W) x 40" (L) x 9" (H)
Retail Box specs	MD20041: 18 lbs MD30002: 27.5 lbs
Warranty	" The M+D Limited Warranty covers your M+D Crutch against manufacturing defects for one year from the date you received your product. Warranty doesnt cover damage caused by accidents or unauthorized modifications."
Weight Capacity	300 lbs