

Capella



Future Mobility Healthcare Future Mobility Healthcare Future Mobility Healthcare



FUTURE

MOBILITY PRODUCTS

Toll-free: 1-855-938-8873 (1-855-9-FUTURE), Fax: 716-783-9236, www.futuremobility.com




The Capella is the most innovative dynamic tilt-in-space wheelchair on the market today. Developed for individuals requiring adjustability to redistribute their body weight for optimal positioning purposes, the Capella offers 20 degrees of tilt, standard on-chair adjustable seating system, incredible range of seat-to-floor heights and a creative streamlined design. The Capella allows the user to function independently.





Adjustable Axle plate position

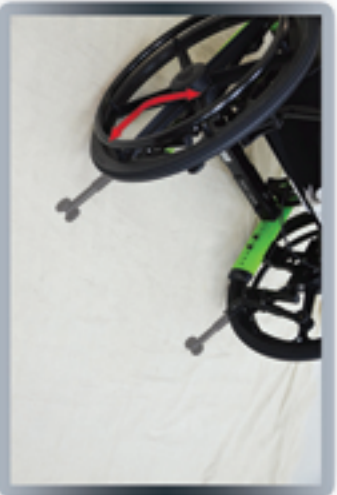

Locking flip up anti tippers



Adjustable front fork angle



Pin release arm allows easy height adjustment of T-style arm



Axle plate allows many different seat to floor heights and center of gravity adjustment
Quick release axle (optional)



Wheel lock available in both push or pull to lock options
Brake extension built into handle



Adjustable back cane angle (83° to 111°)



Chair width adjustable from 16"-20" or 20"-24"



20 degrees of tilt



Capella

Toll-free: 1-855-938-8873 (1-855-9-FUTURE), Fax: 716-783-9236, www.futuremobility.com

SPECIFICATIONS

**MADE IN
CANADA**

Seat Width	16"-24"
Seat Depth	15"-20"
Overall Width	Seat width + 10.75"
Seat to Floor	13"-19"
Rear Wheels	16", 20", 22" or 24" composite urethane or spoked wheels
Caster Size	5", 6", 7", 8" urethane
Back Style	Orion back
Adjustable Back Angle	-7°, 0°, 7°, 14°, 21°
Arm Style	Full length fixed height t-style (gel or vinyl) Full/Desk length adjustable height pin release t-style (gel or vinyl) Full length adjustable height flip back (gel or vinyl)
Front riggings	Swingaway, dual swingaway or dual elevating legrests
Weight Capacity	350 lbs (159 kg)
Product Weight	40 lbs
Frame Material	Lightweight Aluminum
Frame Colours	Black
Lower Frame Colours	



Wet Black Cobalt Blue Purple Grape Green Orange Fire Red Yellow Sky Blue Silver Lining Silver Vein

FEATURES

- 20 degrees of tilt
- 21 degrees of recline (optional)
- Standard seating system consisting of: two way contoured padded fixed headrest, Orion back with 1" depth adjustability and Prism Ideal Cushion
- Aluminum seat pan
- Seat width can increase from 16" to 20" or 20" to 24" as well as seat depth can increase by 1" without the need of changing frame members.
- Multi adjustable axle plate
- Angle adjustable caster fork stem
- Removable rear wheels with stainless steel or plastic coated handrims
- 60 or 70 degree dual swingaway front riggings
- Flip up rear anti tippers
- Padded laptray
- Aluminum or composite footplates
- Heel loops
- Adjustable seat pan (Optional)

WARRANTY

The Capella is warranted to be free from defects in materials and workmanship for a period of one (1) year from the date of purchase and lifetime for the frame. Upholstery for 90 days, under normal usage by original purchaser.

