



FUTURE
MOBILITY

ORION II
by FUTURE MOBILITY



TILT
WHEELCHAIR

THE FUTURE OF MOBILITY – WORKING FOR YOU

"Assembled, Owned and made in North America." In 2000 Future Mobility began designing, manufacturing, and then distributing wheelchairs for the North American Home Healthcare market. Our Founder's mission is to help the world by designing high quality, durable mobility products for the North American marketplace.



ORION II

ORION II 500

ORION III

SUPERKIDS



ORION II

by FUTURE MOBILITY

TILT-IN-SPACE

Future Mobility's Orion II is one of the most innovative, popular tilt wheelchair on the market today. It was developed for active individuals who need to redistribute their body weight to achieve optimal position.



COLOR OPTIONS




A STANDARD WEIGHT LIMIT OF 350 LBS.



STANDARD WEIGHT
LIMIT 350 LBS

NO ADDITIONAL COST



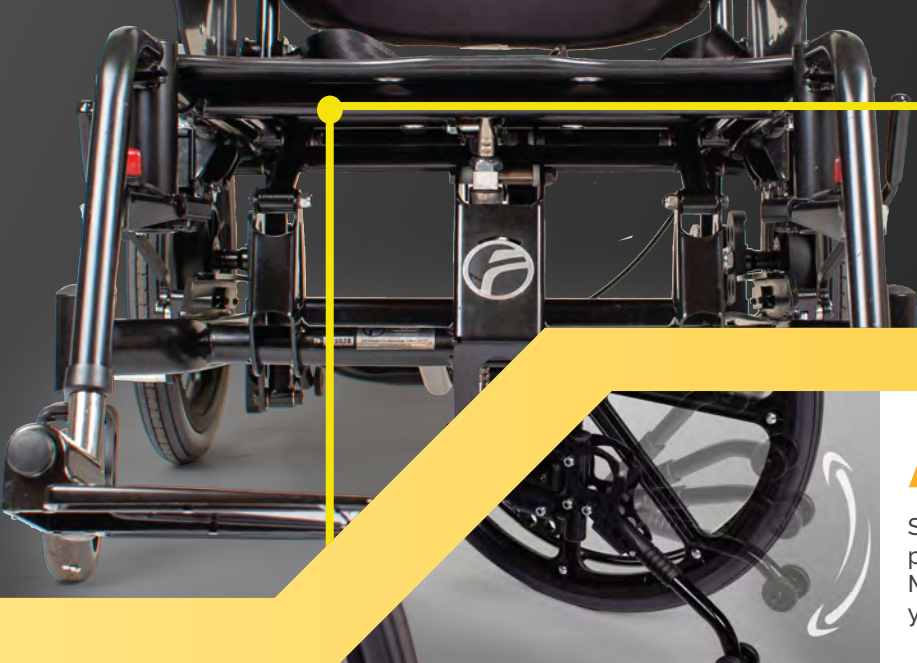


13.25" SEAT-TO-FLOOR HEIGHT

ANTI-TIPPER

Safety is essential. That's why our Anti-Tipper fall prevention system is sold standard with all Future Mobility wheelchairs. The system is designed to keep your loved ones safe from harm.

ARMREST

Future Mobility's wheelchairs are sold with a variety of armrest options. This ensures the end user has the most comfort and support possible when using our wheelchairs. Our Armrests are available in Flip-up, height adjustable, and/or fixed style.

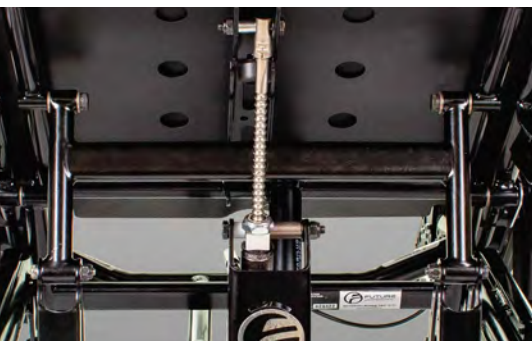

POSI-LOCK

Located underneath the chair, the Posi-Lock allows for a gentle and slow tilt of the end user. It is lockable at different angles.

SWINGAWAY FRONT RIGGINGS


Swingaway Front Riggings give the end user the option of swing-out or dual swing-in and swing-out. Ideal for transporting and seating the end user.



SPECIFICATIONS:



SEAT TO FLOOR
13 1/4" - 20"

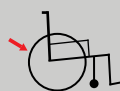

REAR WHEEL
12", 20", 22" OR 24"
COMPOSITE URETHANE



SEAT WIDTH
15" - 24"



SEAT DEPTH
15" - 24"



OVERALL WIDTH
SEAT WIDTH +9.5"



CASTER SIZE
4", 5", 6", 7" AND 8" URETHANE



FRONT RIGGING
PINSTYLE, SWINGAWAY
DUAL SWINGAWAY
OR DUAL ELEVATING LEGRESTS



PRODUCT WEIGHT
65 LBS (29KG)



**ADJUSTABLE
BACK ANGLE**
-7°, 0°, 7°, 14°, 21°



BACK STYLE
ORON BACK



ARM STYLE
ADJ. HEIGHT T-STYLE
ADJ. HEIGHT PIN RELEASE T-STYLE
ADJ. HEIGHT FLIP BACK



WEIGHT CAPACITY
350 LBS (159 KG)

ADP

CODE
WMA5FM005

HCPCS

CODE
E1161

WARRANTY:

The Orion II is warranted to be free from defects in materials and workmanship for a period of one (1) years and the Orion II 500 for period of one (1) year from the date of purchase. The frame has a lifetime warranty and the upholstery for 90 days, upon normal usage by the original purchaser. Power tilt batteries are warranted for six (6) months while the control box, battery charger, actuator and hand control are warranted for one (1) year. If within this warranty period the product shall be proven to be defective, such product shall be repaired or replaced, at Future Mobility Healthcare Inc. discretion.

USA



1 Buffalo River Place, Buffalo,
NY, 14210



1-855-938-8873



www.futuremobility.com



FUTURE MOBILITY

CANADA

6750 Professional Court, Mississauga,
Ontario, L4V 1X6



1-888-737-4011



www.futuremobility.ca

