

Delta™ Ultra-Light 1000, Full-Electric Bed



FIGURE A



FIGURE B

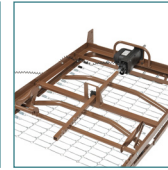


FIGURE C



FIGURE D



INSET E

- The Delta 1000 is truly universal because the bed frame can be used with Drive's and most other manufacturers' bed ends, old and new
- The transition box mounted on the foot section is key to interchangeability (Figure A)
- The transition box allows you to change the rotation of the high/low shaft
- Headboard is taller than footboard to better fit home decor (Figure B)
- The wood grain panels are break and scratch resistant
- Foot section weighs 50% less than conventional foot sections to make deliveries safe and easy
- Once installed, the high/low shaft will stay with the bed to eliminate loss or forgetting on future installations (Figure C)
- Labels and color-coded springs assure accurate installation of side rails for patient safety (Figure D)
- Unique motor is completely self-contained to reduce weight and noise
- The motor assembly can be installed and removed with patient in bed
- Bed, side rails and mattress safety-tested by SGS, a nationally recognized testing facility Satisfies all applicable FDA entrapment guidelines
- The motor is UL approved
- In the event of a power failure, a 9 volt battery, located in the hand control, can lower the head and foot sections nine times No crank needed
- New and improved hand pendant has large, easy-to-use controls (Figure E)
- Bed ships in two cartons
- Channel frame construction provides superior strength and reduced weight
- Zinc-coated spring deck
- Extension kit available (model #15030EXTKIT) extends length from 80" to 84"
- Patent Pending
- Weight Capacity 450 lbs

PRODUCTS	Description
15033	Ultra-Light Full-Electric Bed, 1/cs
15033BV-FR	Ultra-Light Full-Electric Bed with Full-Length Side Rails (item #15001ABV), 1/cs
15033BV-HR	Ultra-Light Full-Electric Bed with Half Rails (item #15201BV), 1/cs
15033BV-PKG	Ultra-Light Full-Electric Bed with Full-Length Side Rails (item #15001ABV) and 80" Innerspring Mattress (item #15006), 1/cs
15033BV-PKG-T	Ultra-Light Full-Electric Bed with Full-Length Side Rails (item #15001ABV) and 80" Therapeutic Support Mattress (item #15019), 1/cs
15033BV-PKG-1	Ultra-Light Full-Electric Bed with Half Rails (item #15201BV) and 80" Innerspring Mattress (item #15006), 1/cs
15033BV-PKG-1-T	Ultra-Light Full-Electric Bed with Half Rails (item #15201BV) and 80" Therapeutic Support Mattress (item #15019), 1/cs
15033BV-PKG-2	Ultra-Light Full-Electric Bed with Full-Length Side Rails (item #15001ABV) and 80" Foam Mattress (item #15007), 1/cs

SPECIFICATIONS	Description
Bed Height Minimum Without Casters	12.75"

SPECIFICATIONS	Description
Bed Height Maximum With Casters	21.5"
Unit Dimensions	36" (W) x 88" (L)
Weight	169 lbs.
Carton Shipping Weight	182 lbs.
Weight Limit	450 lbs Maximum weight limit includes weight of occupant, all accessories such as mattress, side/assist rails, trapeze, any personal affects, etc.
Warranty	(Original Purchaser) 5 Years on Motor and Hand Control, Limited Lifetime on Welds and Frame, 1 Year on All Other Parts and Components
ACCESSORIES	Description
15030HBC	Bedding in a Box - contains 1 pillow case, 1 fitted sheet, 1 cover sheet, and 1 blanket., 1/ea, 5/cs
15030HBL	Fitted Sheet - Contains 2 Fitted Sheets, 2 sheets per bx/ 10 bxs per cs
PARTS	Description
15544H	Head Board
15544F	Foot Board
15544H-F	Set of both
15541	Head Spring Section
15549	Foot Spring Section
15544HPF	Head Board Panel Frame, 1/ea
15544FPF	Foot Board Panel Frame, 1/ea
15544HGB	Head Board Gear Assembly, 1/ea
15544FGB	Foot Board Gear Assembly (w/o crank), 1/ea
15005P-1040	Sleep surface hook, 1/bx
15033MBN	Motor Bracket (Before serial # 1S120329900...), 1/ea
15033MBN-DELTA	(After serial # 1S120329900...), 1/ea
15544GBMB	(After serial # 1S120329900...), 1/ea
15033MHLN-MGBA	Motor Hi/Low & Gear box converter * Must order parts 7b & 7c (After serial # 1S120329900...), 1/ea
15033MHLN-GB	Gear Box Converter Alone * Must order parts 7b & 7c, 1/ea
15033MHLN-S	Mounting Screws, 1/bag
15033MHLNC-N	Motor Shaft (After serial # 1S120329900...), 1/ea
15044CSET	Set of 4 Casters for use with Competitor™ Bed (15560), (2 of ea)
15044CNL	Caster without lock, 1/ea
15044CL	Caster with lock, 1/ea
15030P-1001	Caster Receptacle, 1/ea
15030-FRLINK	Frame Linkage, 1/ea
15033HC-G	Hand Control, 1/ea
15030PN-1035	Head Actuating Bar, 1/ea
15030PN-1036	Foot Actuating Bar, 1/ea
15044P	Foot Board Panel, 1/ea
15544HP-1	Head Board Panel, 1/ea
1503UBS	Hi/Low Shaft
1503C	Collar, 1/ea
15037MN	Motor, 1/ea
15030P-1006	Motor Mount (Okin), 1/bx
15030-7RMM	For Serial No. 7R, 9R use, 1/ea
15030MC	Motor Cap, 1/ea
15030-MS	Motor Strap, 2/pkg
15005P1018	Bed Locking Pin, 1/ea
15544C	Hand Crank, 1/ea
15038MO	Motor, 1/ea
15030UBS	Hi-Lo Drive Shaft (with Collar), 1/ea