


# Etac Molift Quick Raiser 2

Sit-to-stand passive lift with electrically spreadable leg base.

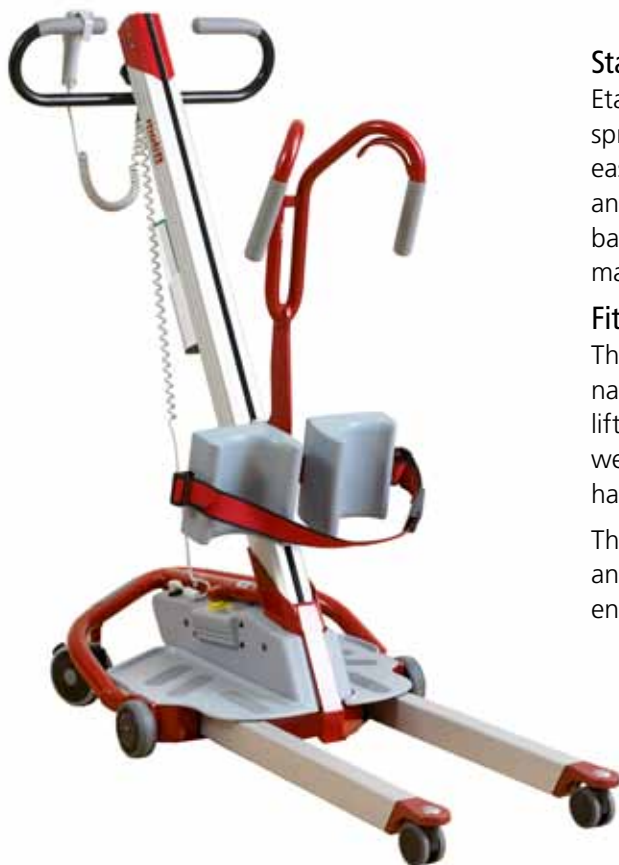

## Stable and easy to maneuver

Etac Molift Quick Raiser 2 with electrically spreadable leg base is compact, stable and easy to maneuver. The unique wheel position and near perfect weight distribution of the base provides a good turning radius and easy maneuverability even in narrow spaces.

## Fits almost everywhere

The inclined straight column imitates the natural pattern of movement and fulfills the lifting needs of users with some degree of weight bearing capacity. Molift Quick Raiser 2 has a safe working load of 160 kg (352 lbs)

The low base, 10 cm (4") fits under most beds and chairs whilst the electrically spreadable legs enables accessibility almost everywhere.

*Find detailed information at page 12–13.*

## Features



### Electrically spreadable legs

The perfect alternative when the lifter is used together with chairs with no accessibility underneath, such as electric wheelchairs and recliners.



### Natural pattern of movement

Inclined straight column imitates a natural pattern of movement.



### Multiple adjustable knee pads

The knee pads are height adjustable & can be turned around for better fitting & to give alternative support for the legs.



### Integrated service software

Automatic service notification system determines when the lifter requires service and maintenance.

## Features – Molift Quick Raiser



### Environmentally friendly NiMH battery

All Molift lifts are equipped with an environmentally friendly NiMH-battery.



### Battery indicator

All Molift lifts have an integrated battery indicator that provides audio and visual notification.



### Unique wheel positioning

The lightweight construction combined with the unique wheel positioning, provides exceptional maneuverability.



### Safety features

Numerous safety features are standard, including emergency stop, overload protection, electrical and mechanical emergency lowering.



### Power-save function

The electronics has a power-save function which makes the system go into sleep-mode after a period with no activity.

## Accessories



### Molift Quick Raiser adjustable footboard

Footboard for adjusting the height of Molift Quick Raiser foot base, ideal for users of smaller height or children. Has a SWL of 75kg (165 lbs)  
Item no: 0451000



### Molift Quick Raiser Sandals

Designed to ensure the feet are held in place during lifting.  
Sandals Small: 0452100  
Sandals Medium: 0452200  
Sandals Large: 0452300



### Extra soft pads for knee support

For patients with extra sensitive legs  
Soft pads Small: 0455111  
Soft pads Large: 0455114

## Options



### Molift Quick Raiser Arm V-Shape

Is the original lifting arm. It is easy to use and fits most users.  
Item no 0458107



### Molift Quick Raiser Arm U-Shape

The newly developed lifting arm has more space between the attachment points for the sling.  
Item no 0458010



### Molift Quick Raiser Knee Support Small

Item no 0455000



### Molift Quick Raiser Knee Support Large

Item no 0455001