

# NON-POWERED Mattress Overlays

Linked by technology, ROHO Mattress Overlays provide unique patient care from prevention to treatment.



## ROHO® DRY FLOATATION® Mattress Overlay System

Create the ideal environment of protection and healing with the original non-powered, air-filled, adjustable, zoned, reactive mattress overlay.

Don't let its simplicity fool you: The performance and efficiency of this overlay system are unequalled in the industry. Comprised of soft, interconnected air-filled cells, the DRY FLOATATION mattress overlay conforms precisely to the contours of the body, helping soft tissue to maintain its shape in order to maximize the healing process.

**Size:** MATTSYSC: 75" x 33" x 3.25" (190.5cm x 86cm x 8.5cm)

**System weight:** 32 lb. (14.5 kg) for the complete system, or 8 lbs. (3.6 kg) per section

**Weight limit:** None if sized to user.

**Warranty:** 2-year limited warranty




**U.S. Medicare Code:**  
E0371 Group 2 Support Surfaces

## PRODUCT FEATURES

- With 720 air cells, this mattress overlay allows for greater conformity and better pressure distribution. Each cell is specially treated to minimize friction and resulting shear.
- The four sections of the mattress overlay are independently adjusted, providing a low-pressure environment over the entire contact area of the body.
- The DRY FLOATATION system provides prevention and therapeutic coverage for all pressure ulcers including deep tissue injuries. The ability to "float" in the mattress overlay promotes healing and reduces soft tissue deformation.
- Using laser Doppler velocimetry in a yearlong study<sup>1</sup>, it was proven that blood flow increased at least four times the base line value when a patient was immersed in the DRY FLOATATION mattress overlay.
- Drainage holes in the base of the mattress overlay allows for moisture to drain away from the body.
- This overlay can be used in any clinical setting. It also is ideal for home, allowing earlier discharge for some individuals.
- There are no hidden costs. This therapeutic non-powered system is a cost effective mattress overlay.

<sup>1</sup> Jacobs, Mary Ann, MSN RN. "Compression of Capillary Blood Flow Using a Regular Hospital Bed Mattress, ROHO Mattress, and Mediscus." Rehabilitation Nursing 14:5 (September/October 1989).

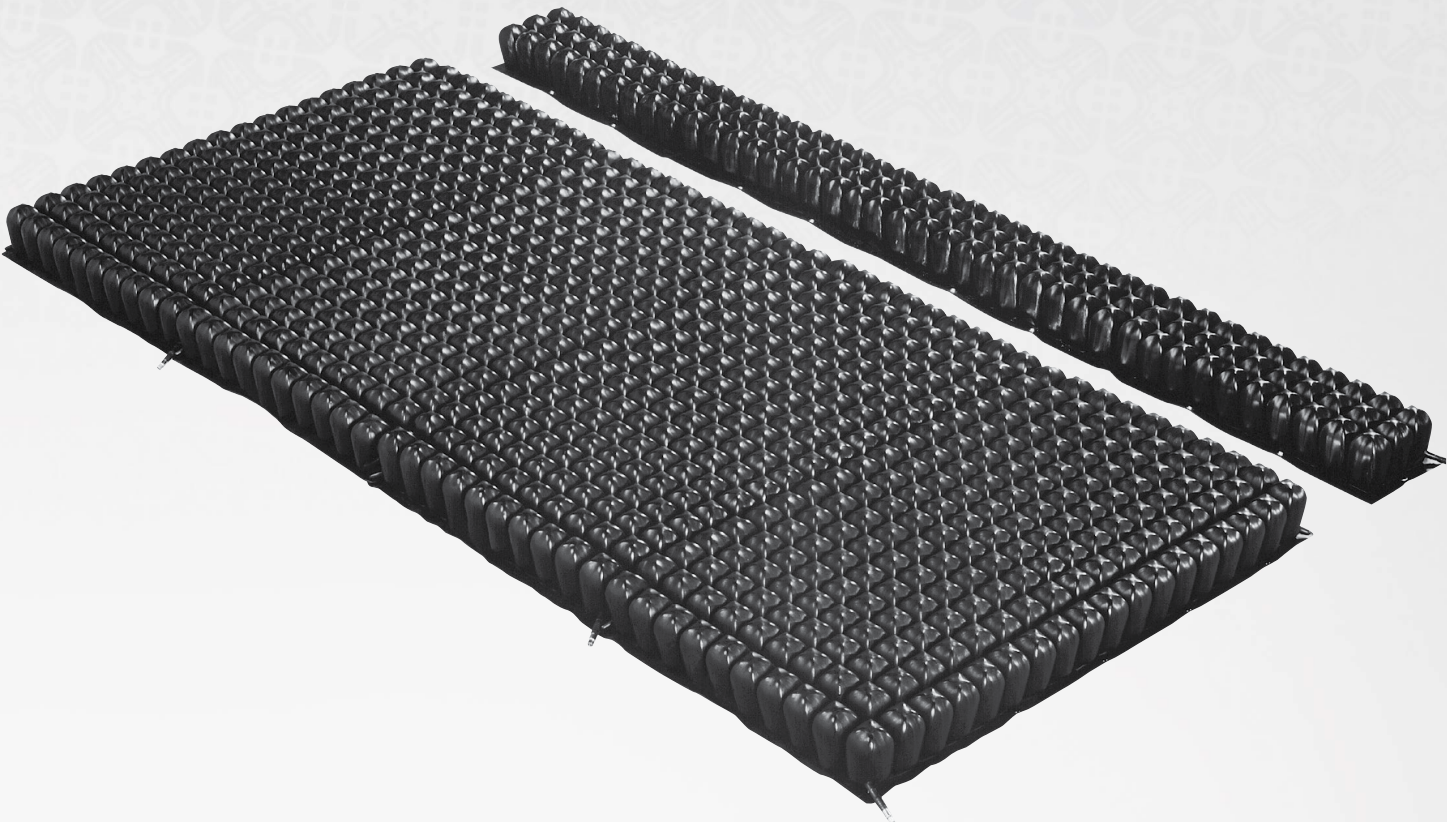

## ROHO® DRY FLOATATION® Bariatric kit

Cell sections extend the width of ROHO DRY FLOATATION Mattress Overlay to fit bariatric bed sizes.

Works in conjunction with a standard size ROHO DRY FLOATATION Mattress Overlay to widen the mattress for use with bariatric bed frames. Add-on sections allow the bariatric individual to experience the therapeutic benefits of the ROHO DRY FLOATATION technology. No user weight limit.

### Size:

BARI39SYC+ 75" x 39.5" x 3.25" / 190.5cm x 100.5cm x 8.5cm  
BARI42SYC+ 75" x 41.5" x 3.25" / 190.5cm x 100.5cm x 8.5cm  
BARI48SYC+ 75" x 48.5" x 3.25" / 190.5cm x 138.5cm x 8.5cm  
BARI54SYC+ 75" x 54.5" x 3.25" / 190.5cm x 100.5cm x 8.5cm  
BARI60SYC+ 75" x 60" x 3.25" / 190.5cm x 100.5cm x 8.5cm

**System weight:** 45 lbs. (20.4kg) based on BARI48SYSC.

**Weight limit:** None if properly sized to the individual.

**Warranty:** 2-year limited warranty



ROHO mattress overlay system is available with your choice of a reusable or enclosure cover.



**U.S. Medicare Code:**  
E0371 Group 2 Support Surfaces