

PERFECT FOR PATIENTS
PERFECT FOR CAREGIVERS





CREATING INNOVATION AND DIFFERENCE

CoZy

The CoZy features 4 zones for catering wounds with capability for automatically customizing mattress pressure to suit individual patient's physique by patient's height and weight.

Cozy Mattress is a medium to high risk system for patients who are at risk of developing Pressure Ulcers according to your sound clinical judgment.

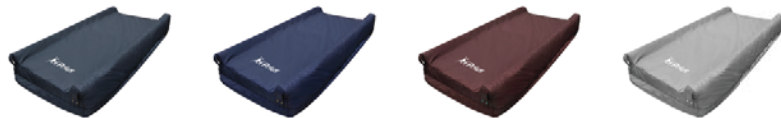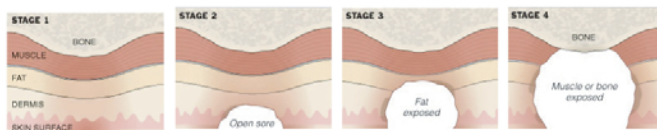

Various top covers are designed and available to satisfy your individual requirement in terms of therapy, prevention and patient comfort.

Mattress Control Unit

- Dimension: (L)12.3"x (W)6.7"x (H)4.9"
- Weight: 14.9 lbs
- Rated Voltage: 110-120 volts
- Max Current: 5A
- Rated Frequency: 60Hz
- Min/Max Pressure: 15-40 mmHg +/- 5 mmHg
- Static / Pulsation Function
- Protection Type: Class I, Type BF
- 10 Digital Adjustable Comfort Control
- Bottom Out Alarm: Visual / Audio
- Leaving Bed Alarm: Visual / Audio
- Power Failure Alarm: Visual / Audio
- User-Friendly LCD Control Panel
- Fowler Function
- Auto Function
- Multi-zone Therapy
- Pulsation

Mattress Replacement

- Model No: CoZy
- Dimension: (L)80"x (W)36"x (H)8"
- Weight: 32.6 lbs
- Weight Capacity: 600 lbs
- Cell Material: Blue Nylon with PU Coating
- Cell Number: 4 cells
- Cover Material: Blue Nylon with PU coating finish
- Base Material: Black Polyester with PVC Coating
- CPR Quick Release System

Toll Free:
877-367-5507
TAMPA, FL OFFICE
TEL: 813-854-5464



CoZy



CoZy

provides you, the carer, with peace of mind
and your patient, with world class care & comfort


Perfect for
patients
and
caregivers

Low friction
base cover
provides
easy contour
with bedframe

Low Air Loss
manages
moisture
and
temperature



Incorporates with
automatic sensing
technology for
bottoming out
prevention and
pressure
adjustment

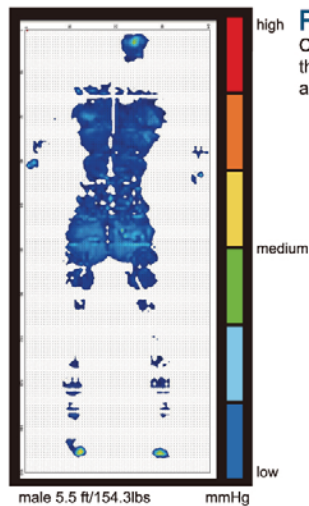

PRESSURE MAPPING

CoZy provides constant low pressure by redistributing the patient's pressure over a greater area than a standard mattress.



- The CoZy mattress replacement can accommodate patient's weight up to 600lbs
- Automatic Fowler boost for upright positioning without loss of therapeutic benefit
- Quilted top cover is highly vapor permeable and oversized providing a low friction and low shear against complications due to moisture
- Convoluted foam base provides for additional safety for patient
- Modular design ensures ease of use and cleaning
- Microprocessor controlled for simplicity, reliability and ease of use
- Reduces risk of patient entrapment and falls



The Low-Pressure Alarm is triggered when the air pressure of the mattress is lower than the default setup.

Each pressure zone can be adjusted independently for customized therapy.

10 adjustable comfort settings provide different levels of comfort depending on the patient's therapy requirement.

The POWER FAILURE LED is triggered during power failure, providing a safer environment for patients.

Auto profile sets mattress pressure based on patient's height and weight.

Capability of customizing individual zone pressure.

CoZy
is your partner
in protection
and treatment
of pressure ulcer



Cable management
keeps
electrical cables
off the floor



CPR quick
release system for
rapid deflation



Hanging bracket
fits on
all type of beds