


Q700M




QUICKIE®

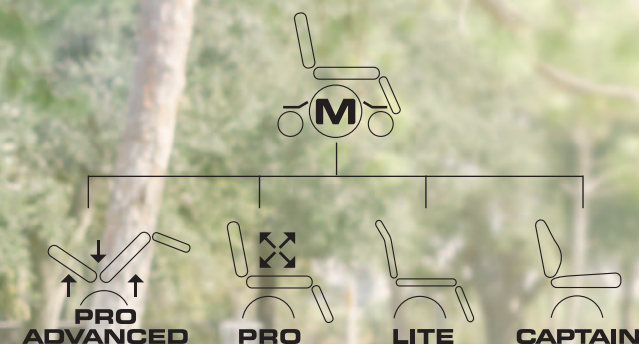

Q700M

THE ULTIMATE
EXPERIENCE

Q700M

Designed to meet your demanding needs.
The Q stands for quality.

ULTRA-RELIABLE MOTORS

HIGH TORQUE, 4-POLE
POLISHED METAL GEARS
SPEEDS UP TO 8 MPH

ULTRA LOW 16.5" STFH

OCCUPIED TRANSIT TESTED TO ISO-7176

SOLID, SINGLE SHEET METAL CHASSIS

CLIMB CURBS UP TO 3"

SIX WHEEL DRIVE SUSPENSION

RELIABLE

Every Q700 M is built to keep you moving. The chassis is formed from a single sheet of metal, and its steel components receive automotive undercoat plating. The industry's best motors undergo additional end-of-line testing before leaving our factory.

MANEUVERABLE

True mid-wheel drive allows the Q700 to turn on its own axis, so you enjoy intuitive driving and excellent maneuverability in tight, compact spaces.

PERFORMANCE

With its independent movement of all six wheels, the Q700 M suspension gives you the ability to safely crawl up and down steep obstacles, keeping all six wheels on the ground.

ELEVATE AND DRIVE WITH C-ME®

C-Me enables you to elevate up to 12" while traveling up to 3 mph! You can lift or descend at the touch of a button with no need to stop while transitioning. Our suspension doesn't lock out when you drive with C-Me, giving you a smooth and comfortable ride from start to finish.





SEDEO® PRO ADVANCED™

Designed for comfort and aesthetics, the SEDEO PRO ADVANCED covers the most demanding clinical needs. With Biometric Repositioning and a full range of programmable seating positions, it is the best ergonomic seating from QUICKIE.



SEDEO® LITE™

The SEDEO LITE is an adjustable seat frame with a sleek look. Its minimal design alleviates weight capacity concerns when using aftermarket equipment. Traditional back canes are width-adjustable from 15" to 20" and easily fit JAY® backs. An array of options cover mild to moderate clinical needs, and intuitive adjustability meets individual requirements.



SEDEO® PRO™

Truly modular in design, the SEDEO PRO provides a new level of functionality, comfort and design. For moderate to complex clinical needs, the PRO combines a highly adjustable seat frame with an excellent choice of postural support options.



Captain's Seat

The Captain's Seat provides a relaxed sitting position and features a semi-recline, fold-down back, adjustable seat depth, and convertible seat pan that allows you to use your preferred cushion. Available with a 50° tilt and power centermount leg support.

Q700 M SEAT FRAME COMPARISON


	SEDEO PRO ADVANCED	SEDEO PRO	SEDEO LITE	CAPTAIN SEAT
Cryo™ Fluid	Available on JAY Syncra Seating	Available on JAY Syncra Seating	Available on JAY Fusion	Available on JAY Fusion
Seat Width	16" to 22"	16" to 26"	15" to 20"	16", 18", 20"
Seat Depth	17" to 22"	14" to 22"	16" to 22"	16" to 22"
Back Support Height	21" to 30"	21" to 30"	18" to 30"	25"
Biometric Repositioning	•			
Memory Seating	•			
Powered Recline	165°	150°		
Anti-shear Back Support	7"	4" (optional)		
Pressure Relief Alert Reminders (PRAR)	•			
Tilt Angle	50°	50°	50°	50°
Elevate with C-Me	12"	12"	12"	
JAY & WHITMYER	•	•	•	•
Remote Seating App	•			

Q700 M See order form as specifications may vary based upon configuration			
HPC codes	No Power Seating Options	Single Power Seating Options	Multi Power Seating Options
Group 3 Standard Duty	K0848/K0849	K0856/K0857	K0861
Group 3 Heavy Duty	K0850	K0858	K0862
Group 4 Standard Duty	K0868/K0869	K0877/K0878	K0884
Max. user weight	300 lbs* (450 lbs for HD)		
Seat-to-floor height	16.5" to 19"		
Total width	25"		
Total length	36.5"		
Turning radius	19.8" (Excludes footplates)		
Max. speed	6 mph (Group 3 Standard)	5 mph (Group 3 HD)	8 mph (optional)
Battery size	Group 34 or Group 24		
Overall weight	Starting at 294 lbs	Starting at 318 lbs	Starting at 366 lbs
Drive wheel position	Mid		
Ground clearance	3"		
Wheel options	14" drive wheels; 6" casters		
Max. range	Up to 18 miles (w/Group 24 batteries and 6 mph motors)		
Transit Brackets	ISO 7176-19:2008 Transit Approved		
Warranty	Lifetime Limited		


*C-Me at 3 mph has a max. user weight of 250 lbs.

COLOR OPTIONS

SHROUD COLORS

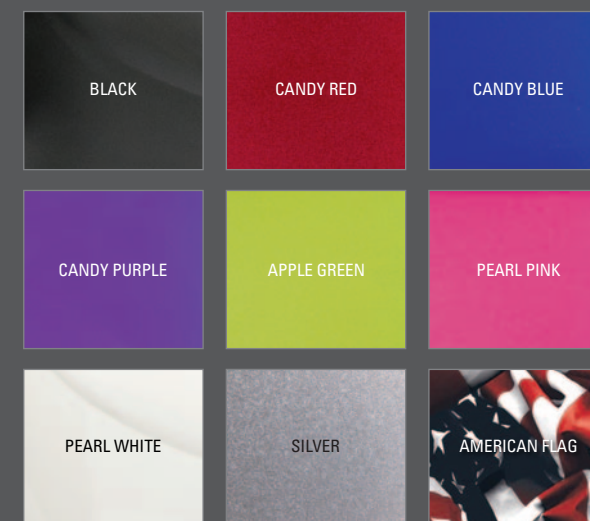

COLOR ACCENTS

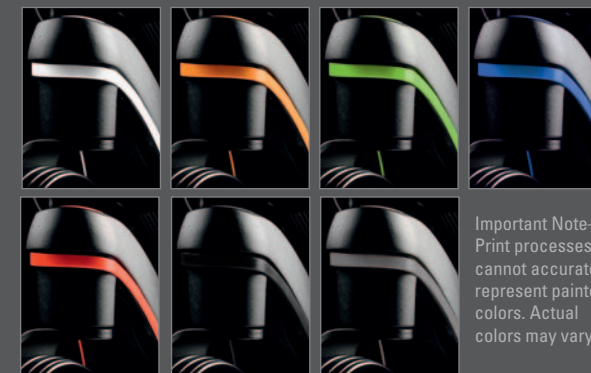

Important Note—
Print processes
cannot accurately
represent painted
colors. Actual
colors may vary.



Web Page



Order Form



Sunrise Medical (US) LLC, Fresno, CA 93727
MK-129475 Rev.H ©03.2022

800.333.4000 www.SunriseMedical.com

