


Salute RDX

The Salute RDX Mattress system provides alternating pressure therapy with micro low air loss to aid patients who are at risk of developing pressure ulcers or have existing stage 1, 2, 3, or 4 ulcers. Caregiver friendly, the Salute RDX is an affordable, therapeutic surface solution for any clinical care setting.

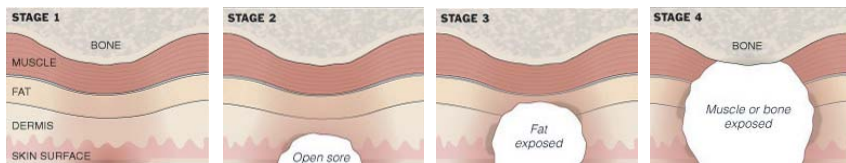

Various top covers are designed and available to satisfy your individual requirement in terms of therapy, prevention and patient comfort.

Master Control Unit

- Dimension: (L)13.5"x (W)7.3"x (H)8.3"
- Weight: 7.1 lbs
- Rating Voltage: 110-120 volts
- Max Current: 0.2A max.
- Rating Frequency: 60Hz
- Cycle Time: 5,10,15,20 mins
- Static Function
- 3 in 1 alternation
- Min/Max Pressure: 12~47 mmHg +/- 6 mmHg
- Air flow rate: 6 liters per minute
- Protection Type: Class I, Type BF
- Low Pressure Alarm: Visual / Audio
- Power Failure Alarm: Visual / Audio
- Maintenance service LED indicator

Mattress Replacement

- Model No: Salute RDX
- Dimension: (L)80"x (W)36"x (H)10"
- Weight: 25.3 lbs
- Weight Capacity: 500 lbs
- Cell Material: Blue Nylon with PU Coating
- Cell in cell design to stay inflated during power outage
- Cell Number: 18 cells
- Cover Material: Blue Nylon with PU coating finish
- Base Material: Black Polyester with PVC Coating
- CPR Quick Release System
- Cable Management System on bottom cover

Toll Free:
877-367-5507

TAMPA, FL OFFICE
TEL: 813-854-5464



Prius

CREATING INNOVATION AND DIFFERENCE



Salute RDX


HAPPYHEEL™ Heel Relief
Patent Pending



Salute RDX

provides you, the caregiver, with peace of mind and your patient, with world class care & comfort

Perfect for patients and caregivers

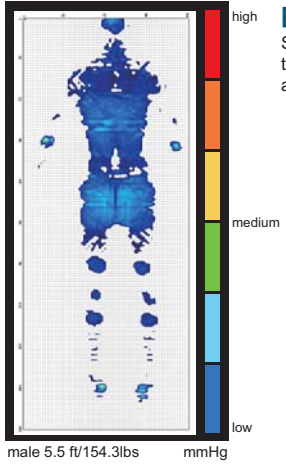

PRESSURE MAPPING
 Salute RDX provides constant low pressure by redistributing the patient's profile over a greater area than a standard mattress.



- Interchangeable air cell design provides easy cleaning and minimizes the proliferation of bacteria
- The Salute RDX employs 8" air bladders with cell in cell design (4.5" of therapy, 3.5" of static air), to provide therapeutic alternating pressure and also allow the mattress to stay inflated in a power outage.
- The Salute RDX replacement mattress can accommodate patient weights up to 500 lbs

HAPPYHEEL™ provides extra softness in the last 5 air cells to provide protection for patient's fragile heels

- Modular design ensures ease of use and cleaning

Nylon covered bottom cover for friction and noise reduction when raising bed's head pan


Effective micro low air loss for moisture management



Slimline pump design is unobstructive

8 adjustable comfort settings provide different levels of comfort depending on the patients therapy requirement.

HAPPYHEEL™ provides further therapy and comfort by decreasing pressure in patient's heel area.

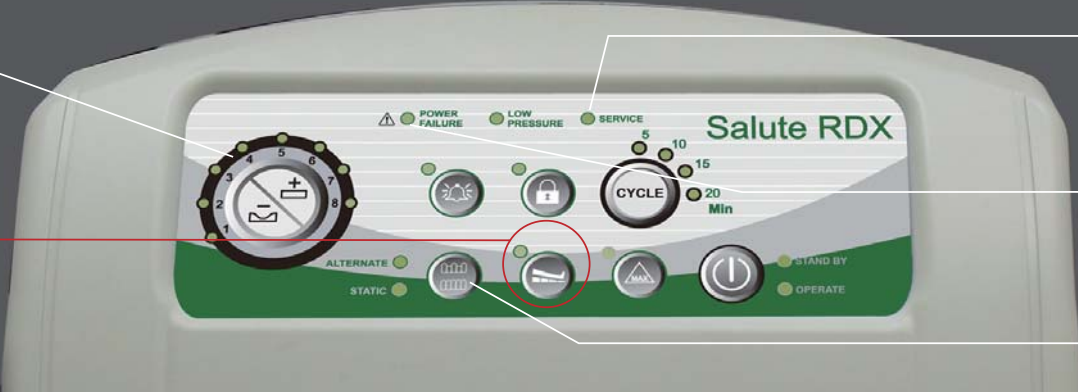

Service LED and alarm indicator is a friendly reminder for advising caregiver to send back control unit to certified technician for annual inspection.

The Power Failure alarm warns in case of power outage.

3-in-1 alternating pressure covers 66% of patient's body (vs. 50% on 2 in 1 systems) to provide a more comfortable yet still effective therapy

Salute RDX is your partner in prevention and treatment of pressure ulcers

Cable management system keeps electrical cables off floor for caregiver and patient safety

CPR quick release system deflates the system within 20 seconds



Spring loaded hooks snap to fit your bed's footboard