

Expanding the Dimensions of Acute Care

Ai1 All in One expands what is possible in a cost effective med-surg bed.

Engineered for Versatility and Expanded Use

Designed for Patient Safety and Security

Enhanced Falls Management

Unique Platform Design

Designed for an expanded experience.

Expanded width

Expanded capacity

Expanded care



Designed for Patient Safety and Security

Safely care for med-surg patients with nurse communications monitoring, all the essential bed features you expect, with a full line of surfaces proven to prevent and treat pressure injuries.



Nurse Call Communications

Enables connectivity that is intuitive to use, access, and activate. Compatible with most nurse call systems.



In-Bed Scale

Accurately monitors and records patient weight and weight history



Safe Siderails

Minimized risk of patient entrapment and pinch points with siderails designed for full compliance with the International Electrotechnical Commission (IEC), the IEC60601-2-52:2009



Convenient, Flat Controls

Convenient access to flat control panels for both patient and caregiver, supporting easy and effective infection control



Safety Stop Technology

Automatically stops bed lowering if contact detected

Enhanced Falls Management

Facilitate purposeful early mobility with an ultra-low bed height and 3-zone bed exit connectivity designed to create a safer environment.



Low Bed Height

Enhance safe ambulation and minimize falls with a lower bed height that fits individual patient needs


Easy-View, Split Siderails

Ergonomically designed to optimize patient visibility while assisting active patient repositioning with easy-grip hand holds



3-Zone Bed Exit

Enhance response time with connected monitoring and alarms

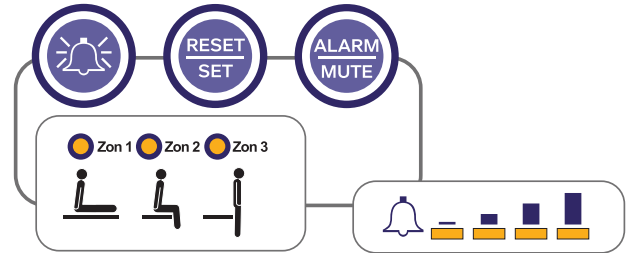

Under Bed Light

Softly illuminates bed silhouette for improved night-time visibility yet minimizes sleep disruption



Unique Platform Design

Expandable bed platform designed specifically to accommodate a wide range of patient profiles. No longer will hospitals be at risk when larger patients are admitted and sent to a med-surg unit. The Ai1 has you covered.


Expanding Care Features

Low bed height

Expandable width

Expandable length

600 lb. /550 lb.
(SWL/Max Pt.)

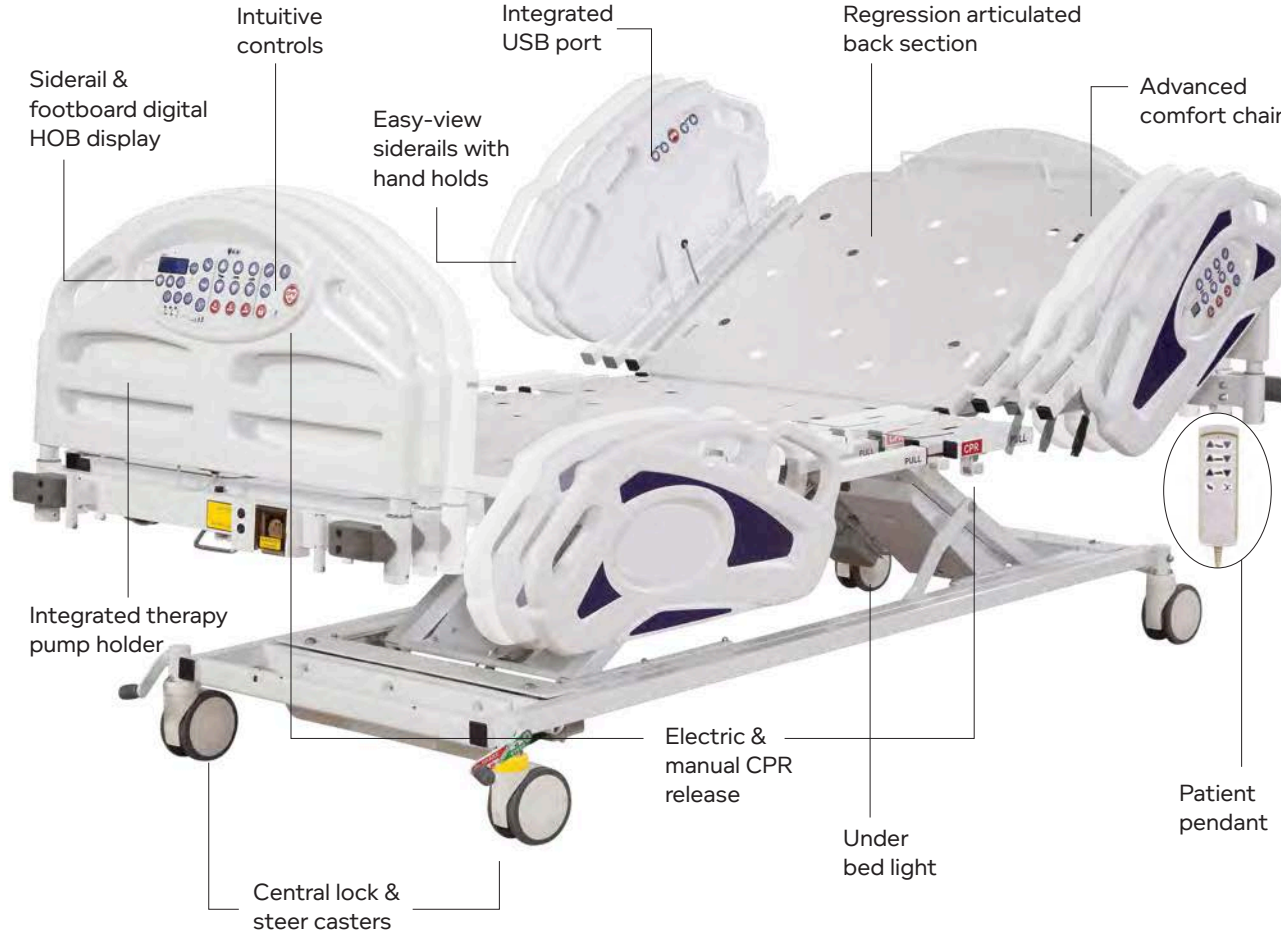

Premium Options

In-bed scale

Nurse call/
communications

3-zone
bed exit alarm

On-board battery

Integrated hospital-
grade outlet

In-rail controls

Safety stop
technology



Bed Specifications

Expandable length	80" - 84" - 88"
Expandable width	35" - 42" - 48"
Safe working load/max. patient weight	600 lbs./550 lbs.
Weight	360 lbs
Height	<10" - 31"
Head of bed angle	0 - 70°
Knee angle	0 - 35°
Trendelenburg/reverse Trendelenburg	14°
Caster diameter	5"
Mattress size	36" - 48" wide 80" - 88" long
Power requirements	115 VAC, 60 Hz, 4 AMP
Safety standard	IEC60601-2-52