



Model Ci



# The New Era of Personal Mobility Has Finally Arrived.

The WHILL Model Ci is a revolutionary approach to power mobility that uniquely boosts user confidence and evokes an active, can-do lifestyle. Its premium performance and superb engineering provides versatility and reliability to experience all-day adventures.

From socializing with friends and family to exploring the great outdoors, Model Ci empowers users to enjoy their world freely both indoors and out.

Powerful and agile, the Model Ci is exhilarating and fun to drive. Its distinctive design allows you to be you.


Pursue life's adventures with the Model Ci Personal EV.






## Refuse to Accept Limitations


With its compact frame, patented front omni-wheels, and responsive mouse controller, Model Ci has the maneuverability to effortlessly navigate confined spaces.

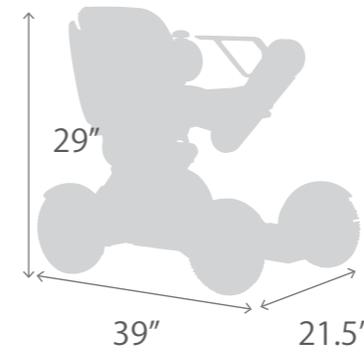


Smaller footprint than traditional power mobility devices.

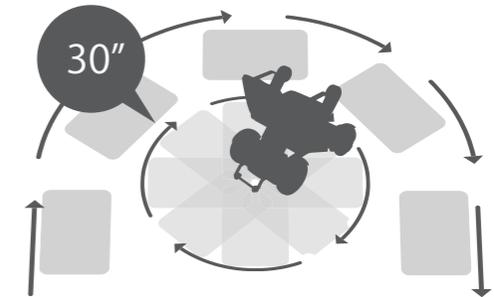

Up to a 45% tighter turning radius than tiller operated mobility devices.

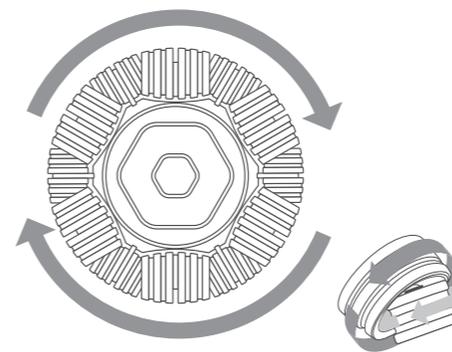

Patented front omni-wheels for tight turning and superb terrain coverage.

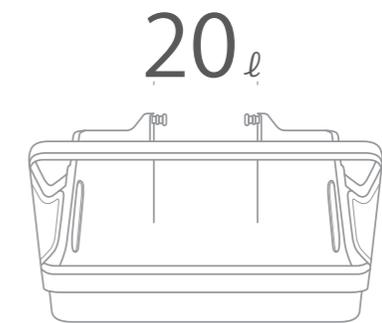

Large under-seat basket to store your belongings and groceries.

# Get Out and Adventure

Equipped with two powerful motors and engineered for superb stability, effortlessly go from indoor activities to navigating trails and surrounding yourself with nature.






## Live Unpredictably

Model Ci quickly disassembles into three lightweight parts for transport in any vehicle with room to spare. Assembly is just as easy to get you quickly moving again at your destination.

Quickly Disassemble

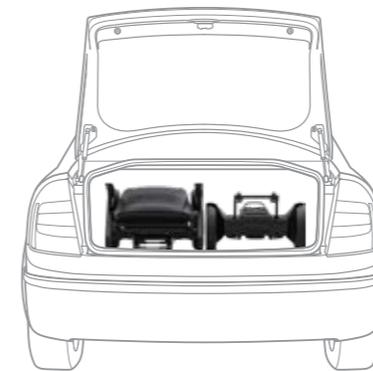


Load Simply



Store Neatly

# Connect Intelligently

Equipped with intelligent innovations, such as the very latest in app controls and mobile data network connectivity, for greater independence and premium support.

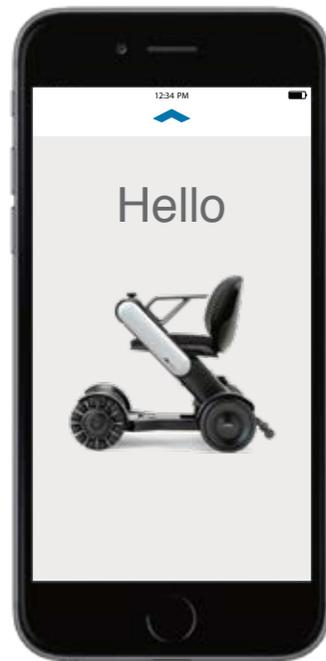

Control remotely and manage speed settings with the iPhone app.



Lock/Unlock your device with the WHILL Smart Key or iPhone app.




Charge your mobile devices with a built-in USB port.





## Drive with Peace of Mind

Keep engaged in your activities and embrace your surroundings by knowing Model Ci will provide the necessary response and feedback to keep you riding safe.

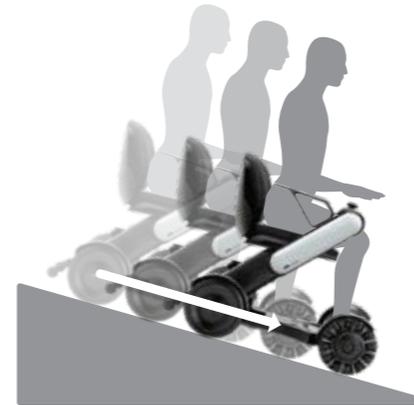

Electromagnetic brakes provide safe stopping, even on inclines.



An anti-sway system keeps you driving straight on side slopes.



A built-in horn notifies others of your presence to avoid accidents.

Audible messages notify you of operational errors.



Arms rotate up for easy access to the seat.

# Customization

Match your personality and style with six color options.



# WHILL Personal EVs

Find the personal EV that best fits your lifestyle at [whill.us](http://whill.us).



Model Ci

Ultra-Portable



Model A

Indoor/Outdoor Mobility  
with AWD



Model M

Next-Gen Power Wheelchair  
with AWD

# Specifications

## Performance Specifications

Driving Range	10 Miles	Turning Radius	29.9"
Maximum Speed	5 MPH	Obstacle Clearance	2"
Drive System	Two Motor Drive System	Ground Clearance	2.25"
Incline Capability	10°	Braking System	Electromagnetic
Controller	WHILL Control System	Max Weight Capacity	250 LB

## Dimensional Specifications

Chair Weight	115 LB	Front Wheel Diameter	9.8" Omni-wheel
Chair Width	16" Seat = 21.5" / 18" Seat = 23.5"	Rear Wheel Diameter	10.5"
Chair Length	38.8"	Rear Tires	Non-Puncture

## Other Specifications

Charging Time	5 Hours	Operating Conditions	5°F to 104°F
Battery Charger	Lithium-ion (25.2V 10Ah)	Tail Lamps	LED Light (Red)
Seat Cushions	16"x16", 16"x18", 18"x18"	Horn	Built-In
USB Port	Single	Lap Belt	Optional

# Accessories to Keep You Moving



Smart Key



Side Bag



Lap Belt



Spare Battery



Spare Charger



Walking Stick Holder

NOTE: Accessories shown are available for Model Ci only.

# Intelligent Personal EVs

[www.whill.us](http://www.whill.us)

WHILL Inc.  
285 Old County Road, #6,  
San Carlos, CA 94070

The WHILL Model Ci is not considered to be a medical device and has not been submitted to the Food and Drug Administration for review or clearance.

740-02270-MKT-01 12/17