



Toll-free: 1-855-938-8873 (1-855-9-FUTURE), Fax: 716-783-9236, www.futuremobility.com

The Super Chair for Kids by Orion is a low maintenance, high performance tilting wheelchair. The Super Chair for Kids by Orion comes with many standard features, including flip up footplates which can be adjusted to an angle comfortable for your clients needs. Our selection of bright colours and a chair developed for flexibility and comfort, makes this wheelchair the model of choice as it grows with its occupant!







Swinging front leg rests either 70, 80 or 90 Degrees and provide the means to tie down when wheelchair is in a vehicle Front and rear transit hooks come standard with the wheelchair with angle adjustable footplates transfer ites can be flipped up SUPER CHAIR FOR KIDS Comes with vinyl arm pad T - Style arm can be easily adjusted up or down. tilt position and is Posi Lock Mechanism locks the wheelchair in controlled by the tilt handle heights and center of gravity adjustment Axle plate allows many seat to floor up or down Rear anti tippers can be adjusted easily tilt chair Tilt handle allows attendant to both push or pull to lock handle Brake extension built into Wheel lock available in

Removable stroller handle

SPECIFICATIONS

- Seat Width 10"- 16" - Seat Depth 11"- 16"

- Overall Width Seat width + 9" MADE IN
- Seat to Floor 17 1/2" - 20" CANADA

- Rear Wheels 12", 20", 22" composite urethane 4", 5", 7", 8" urethane

- Caster Size 4", 5", 7", 8" uretnane
- Back Style Super Chair Kids Back.

- Back Height 19" – 23" (from top of seat pan) - Adjustable Back Angle -6°, 0°, 6°, 12°, 18°, 24°, 30°

- Arm Style Full/desk length adjustable height t- style (gel or vinyl)

Full/desk length angle adjustable cantilever arm (gel or vinyl)

- Front riggings Swingaway or elevating legrests

- Weight Capacity 150 lbs (68 kg)

- Product Weight 39 lbs

- Frame Material Lightweight Aluminum

FEATURES

- 3-40 degrees of tilt


- Standard seating system consisting of: two way contoured padded fixed headrest, Super Chair Kids Back and Super Chair Kids Cushion
- Seat width can increase from 10" to 13" or 13" to 16" as well as seat depth can increase from 11" to 14" or 13" to 16" without the need of changing frame members.
- Aluminum front and rear seat pan
- Multi adjustable axle plate
- Arm heights: 8-1/2" 12-1/2" (adjustable)
- Removable rear wheels with stainless steel or plastic coated handrims
- 70 degree swingaway upgraded front rigging
- Front and rear transit hooks
- Adjustable rear anti tippers
- Aluminum or composite footplates
- Heel loops

WARRANTY

Future Mobility Products warrants its Super Chair for Kids by Orion Wheelchair to be free from defects in materials and workmanship for one (1) year from the date of purchase. The frame is warranted for the lifetime of the original purchaser/user and the upholstery for 90 days, upon normal usage by original purchaser. If within this warranty period the product shall be proven to be defective, such product shall be repaired or replaced, at Future Mobility Products discretion.

FDA CODE: **K103073**

FRAME COLOURS

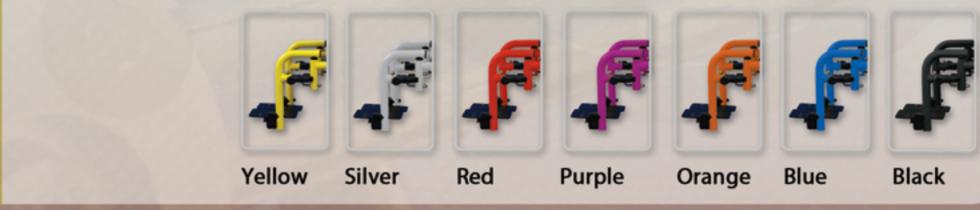

1 Buffalo River Place Buffalo, New York, 14210, USA
Toll-free: 1-855-938-8873 (1-855-9-FUTURE), Fax: 716-783-9236, www.futuremobility.com

