

Multi-Ply 6500 Dynamic Elite Pressure Redistribution Foam Mattress

HCPCS: E0184

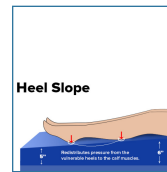
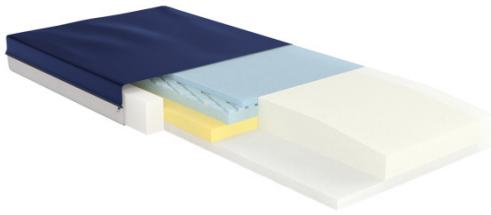


FIGURE A

- High quality therapeutic foam mattress provides pressure redistribution and shear/friction reduction
- Drive's premium therapeutic foam mattress maximizes value and comfort with 4 full layers of latex free foam
- Full width gel infused, visco-elastic foam top layer dissipates heat providing a cool and dry environment to help preserve skin integrity
- Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress
- High density foam perimeters provide solid edge support and facilitate safe patient transfers
- Available with or without raised side rails
- Heel section is sloped with patented LuraQuilt foam offloading pressure from vulnerable heels
- Removable premium, multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Meets flammability standard CFR 1632 and CFR 1633

PRODUCTS	Description
6500-DE-2-FB	6500DynamicElite Mattress 80",1/ea
SPECIFICATIONS	Description
Cover Material	Fluid-Resistant Stretch Top Waterproof Vinyl Bottom
Dedicated Heel Section	Yes
Fire Rating	CFR 16 Part 1633 Federal Fire Code
Mattress Layers of Foam	4
Side Rails	Yes
Warranty	15 Year Limited Foam 3 Year Cover
Weight Capacity	400 lbs
ACCESSORIES	Description
CPU80	Cover,MultiPly/ShearCareSeries,6x36x80, EA