


Choices to help prevent pressure ulcers
and improve outcomes.

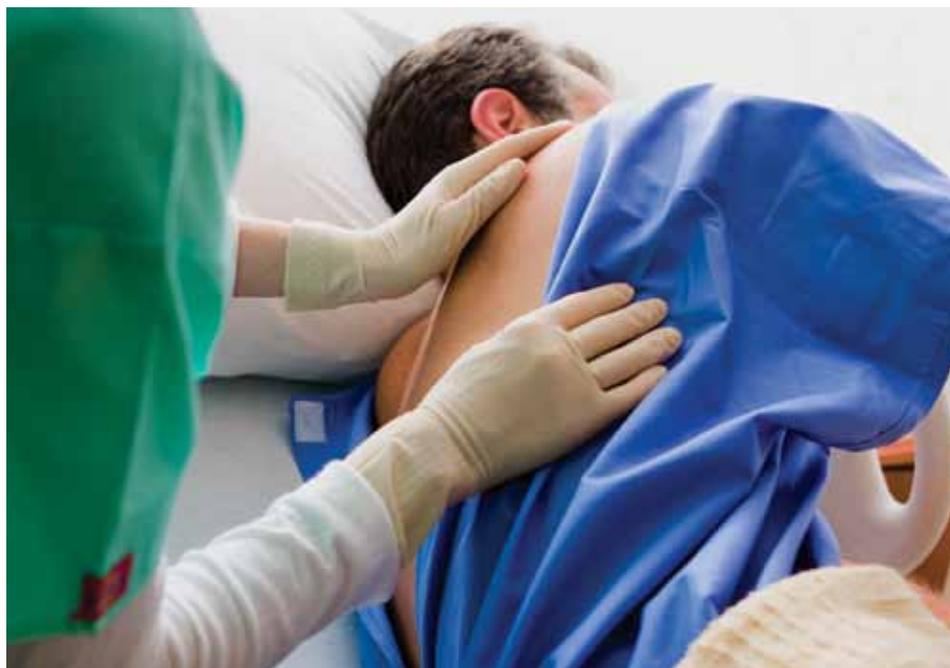

Non-Powered Surfaces Portfolio

Hill-Rom

Enhancing Outcomes for Patients and Their Caregivers.™

Hill-Rom's non-powered surface portfolio, coupled with the Safe Skin® program, can help reduce your pressure ulcer prevalence.

Reactive Non-Powered Surfaces



Hill-Rom® NP50 Prevention Surface

- Single-density hospital grade foam
- Designed to articulate with a frame
- Fluid proof nylon ticking and vinyl base
- Optional bolster/raised firm perimeter
- Accommodates weight up to 350 lbs or 700 lbs on bariatric version

Patient Indications:

- Low risk patients
- Prevention

Hill-Rom® NP100 Prevention Surface

- Multi-zoned foam sections and layers allow for improved pressure redistribution points
- Shear liner for added patient comfort and shear prevention
- Bi-directional stretch ticking reduces friction and noise
- RF Welded seams to prevent fluid ingress
- Accommodates weight up to 350 lbs or 500 lbs on specialty or bariatric versions

Patient Indications:

- Low risk patients
- Prevention

Hill-Rom® Tempur-Pedic® Medical Wound Surface

- Reduces the risk of pressure ulcers by utilizing clinically proven Tempur-Pedic® material
- Makes the patient more comfortable and relaxed within the hospital atmosphere with superior comfort levels known by Tempur-Pedic®
- Accommodates weight up to 500 lbs

Patient Indications:

- Low-moderate risk patients
- Prevention and treatment for Stage I, II

TECHNOLOGY

COMFORT
PREVENTION

TREATMENT



AccuMax Quantum™ VPC Pressure Redistribution Mattress

- Non-powered Controlled Release™ technology
- Multi-zoned and layered with varying ILD's and densities
- Provides uninterrupted pressure redistribution through a series of interconnected air sectors
- No blower box, pump or power required
- Patented Heel Pillow™ with heel slope*
- Accommodates weight up to 500 lbs

Patient Indications:

- Moderate risk patients
- Prevention and treatment for Stage I and II and uncomplicated III's and IV's with multiple turning surfaces

AccuMax Quantum™ Convertible Pressure Redistribution Mattress

- Non-powered Controlled Release™ technology
- Multi-zoned and layered with varying ILD's and Densities
- Provides uninterrupted pressure redistribution through a series of interconnected air sectors
- No blower box, pump or power required
- Optional Alternating Pressure at the touch of a finger
- Easy access airport allows for quick access to control unit for powered therapy
- Patented Heel Pillow™ with heel slope*
- Accommodates weight up to 500 lbs

Patient Indications:

- Moderate risk patients
- Prevention and treatment for Stage I and II and uncomplicated III's and IV's with multiple turning surfaces

Hill-Rom® NP200 Stimulite® Honeycomb Wound Surface

- Multi-zoned, unique aerospace technology, Stimulite® Honeycomb, provides self-adjusting pressure relieving properties
- Maintenance free, stand-alone surface with no power, blower box or additional parts necessary
- Light surface weight of only 24-26 lbs
- Accommodates weight up to 500 lbs

Patient Indications:

- Moderate risk patients
- Prevention and treatment for Stage I and II and uncomplicated III's and IV's with multiple turning surfaces

Frames

Offering a variety of non-powered mattresses designed to fit Hill-Rom® beds and stretchers as well as other flat deck frames.



VersaCare® Bed

The med-surg bed that goes the extra mile. This bed has the lowest height in acute care and offers superior patient comfort and intuitive controls that help caregivers easily operate the system to its fullest capacity.

Non-Powered Surface Compatibility:

- Hill-Rom® NP200 Stimulite® Honeycomb
- AccuMax Quantum™ Convertible/AD
- AccuMax Quantum™ VPC/IS
- Hill-Rom® Tempur-Pedic® Medical Mattress
- Hill-Rom® NP100



CareAssist® Bed/CareAssist® ES Bed

Its contemporary appearance, improved braking system, enhanced patient positioning, and easier-to-use controls, all help make this bed one that meets the needs of everyone in the acute care setting.

Non-Powered Surface Compatibility:

- Hill-Rom® NP200 Stimulite® Honeycomb
- AccuMax Quantum™ Convertible/AD
- AccuMax Quantum™ VPC/AD
- Hill-Rom® Tempur-Pedic® Medical Mattress
- Hill-Rom® NP100



TotalCare® Bed

A versatile bed that benefits patients of all kinds. It gives caregivers the power of touch at the point of care with its clean, simple controls and a large, colorful touchscreen on the siderails.

Non-Powered Surface Compatibility:

- AccuMax Quantum™ Convertible/TC
- Accumax Quantum™ VPC/TC
- Hill-Rom® NP100



Stretchers

Strong, robust stretchers that meet the requirements of the toughest working conditions; designed for mobility, comfort and quality, Hill-Rom stretchers are built for outstanding performance with multiple prevention surface options.

Non-Powered Surface Compatibility:

- Hill-Rom® Tempur-Pedic® Medical Mattress
- Comfortline® mattress
- AccuMax Quantum™ VPC



Flat Deck Beds

Durable flat-deck solutions meeting the flexibility of the med-surg environment.

Solutions for Basic Care™, Advance™, Advanta™ 2 and Century Beds and other legacy frames.

Non-Powered Surface Compatibility:

- Hill-Rom® NP200 Stimulite® Honeycomb
- AccuMax Quantum™ Convertible
- AccuMax Quantum™ VPC
- Hill-Rom® Tempur-Pedic® Medical Mattress
- Hill-Rom® NP100
- Hill-Rom® NP50

Options vary by model.

AccuMax Quantum™ is a trademark of Encompass Group, Inc. Stimulite is a registered trademark of Supracor, Inc.

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

© 2010 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.

Hill-Rom

Enhancing Outcomes for Patients and Their Caregivers.™

USA 800-445-3730
Canada 800-267-2337

www.hill-rom.com