



# ROHO<sup>®</sup> PRODIGY MATTRESS OVERLAY<sup>®</sup> SYSTEM

## SKIN PROTECTION

Protecting an individual's skin/soft tissue is critical for prevention and treatment of pressure ulcers. The ROHO PRODIGY Mattress Overlay System is a non-powered, adjustable and reactive support surface for individuals who require skin/soft tissue protection or have other therapeutic needs. The PRODIGY Mattress Overlay provides individualized pressure, shear and friction management through force redistribution to an individual's shape.

## ADJUSTABLE

Flexible, interconnected air cells can be adjusted to match the individual's body contour, providing customized comfort and skin/soft tissue protection. The PRODIGY Mattress Overlay has three sections for an individualized fit. Each section can be adjusted to address the specific needs of the individual.

## LIGHT WEIGHT

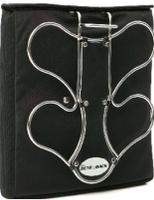
A complete PRODIGY Mattress Overlay System weighs only 9.5 lbs (4.3 kg). The three-section PRODIGY Mattress Overlay is so easy to carry and handle during set up, inflation and adjustment, that it is the ideal choice for protection in any care setting.

## NON-POWERED

No need for concern of power failure or energy consumption costs.



# ROHO® PRODIGY MATTRESS OVERLAY® SYSTEM

APPLICATIONS	FEATURES & BENEFITS	MORE PRODUCTS FROM ROHO
<p><b>THE PRODIGY MATTRESS OVERLAY IS EFFECTIVE FOR INDIVIDUALS:</b></p> <ul style="list-style-type: none"> <li>- who have limited mobility or completely immobile</li> <li>- are at risk of developing pressure ulcers and/or tissue injury</li> <li>- are at risk for skin/soft tissue breakdown</li> <li>- need comfort</li> <li>- need pain management</li> <li>- need a mattress overlay that provides an individualized fit</li> </ul>	<p><b>AIR FLOATATION TECHNOLOGY:</b> Air Floatation technology redistributes forces that can cause peak pressure, shear and friction. The PRODIGY Mattress Overlay is effective for pain management, comfort and individuals who require skin/soft tissue protection and treatment.</p> <p><b>TRANSPORTABLE:</b> The entire PRODIGY Mattress Overlay System weighs only 9.5 lbs (4.3 kg). The air can be released for easy transport between the hospital, care facility, home or travel. No need to buy a new mattress overlay for each setting.</p> <p><b>ADJUSTABLE:</b> Three sections create a multi-zoned support surface for an individualized fit. Each mattress overlay section can be adjusted per body segment to address the specific needs of the individual's shape and size.</p>	 <p><b>ROHO® Hybrid Elite® Cushion</b></p>  <p><b>ROHO® HEAL PAD® Cushion</b></p>
<p><b>SPECIFICATIONS</b></p>		 <p><b>ROHO® MOSAIC® Cushion</b></p>
<p><b>INCLUDES:</b> PRODIGY Mattress Overlay System (3 sections), straps, hand inflation pump, operating instructions and repair kit.</p> <p><b>CONSTRUCTION:</b> Flame-Resistant Polyvinyl</p> <p><b>SIZES:</b> MATTRESS OVERLAY SYSTEM (3 Sections): SYSTEM: 36" x 81.75" x 3.5" (91.5 cm x 207.5 cm x 9 cm) SECTION: 36" x 27.5" x 3.5" (91.5 cm x 70 cm x 9 cm)</p> <p><b>WEIGHT:</b> A COMPLETE PRODIGY SYSTEM: 9.5 lbs (4.3 kg) A Complete PRODIGY System includes three (3) Multi-Zoned Sections</p> <p><b>WEIGHT LIMIT:</b> 300 lbs (136 kg) for complete system</p> <p><b>WARRANTY:</b> 12-month limited warranty</p> <p><b>U.S. MEDICARE PART B CODE:</b> E0197 Approved for Group 1 Support Surfaces</p>	<p><b>FITS MOST BEDS:</b> Can be placed on top of almost any size mattress, including the individual's existing home mattress, hospital bed mattress or care facility bed mattress.</p> <p><b>LIGHTWEIGHT:</b> A complete PRODIGY Mattress Overlay System weighs only 9.5 lbs (4.3 kg). Easy to carry and handle during set up, inflation and adjustment.</p> <p><b>NON-POWERED:</b> Energy-efficient because it requires no power source for use. There is also no noise or heat generated from a pump.</p> <p><b>COVER OPTIONS:</b> Two options designed to fit every ROHO Mattress Overlay System: ROHO® Reusable Mattress Overlay Cover (one piece that tucks in like a fitted sheet) or the ROHO® Mattress Overlay Enclosure Cover (has a top and bottom that zip together to completely enclose the mattress overlay). These two cover options are made in America and are fluid-resistant, flame-resistant, and can be easily removed for cleaning and disinfecting between uses.</p>	 <p><b>Retroback® Back Support System</b></p>  <p><b>ROHO® QUADTRO SELECT® HIGH PROFILE® Cushion</b></p> <p><b>FOR ADDITIONAL INFORMATION, PLEASE CONTACT:</b></p> <div style="border: 1px dashed black; height: 100px; width: 100%;"></div>



100 North Florida Avenue • Belleville, IL 62221-5429 • USA • www.therohogroup.com  
 U.S.: 1 800 851 3449 • 1 618 277 9173 • fax 1 888 551 3449  
 Outside the U.S.: 1 618 277 9150 • fax 1 618 277 6518

©2008, 2010, 2011 ROHO, Inc. The following are registered trademarks of ROHO, Inc.: DRY FLOATATION®, HEAL PAD®, HIGH PROFILE®, Hybrid Elite®, MOSAIC®, PRODIGY Mattress Overlay®, QUADTRO SELECT®, Retroback®, ROHO® and shape fitting technology®.

ROHO, Inc. has a policy of continual product improvement and reserves the right to amend this document. The current version of this document is available on our website at www.therohogroup.com.  
 LIT002120, Rev. 08/11, SP 08/11



MADE IN USA  
 PRODIGY Mattress Overlay

